

# openheart A new electrode placement method for obtaining 12-lead ECGs

Gabriel M Khan

**To cite:** Khan GM. A new electrode placement method for obtaining 12-lead ECGs. *Open Heart* 2015;2:e000226. doi:10.1136/openhrt-2014-000226

► Additional material is available. To view please visit the journal (<http://dx.doi.org/10.1136/openhrt-2014-000226>)

Received 2 December 2014  
Revised 17 April 2015  
Accepted 6 May 2015



CrossMark

Dr Gabriel M Khan MB, BCh, MD (Queen's Belfast), FRCP (London) FRCP C, FACP, FACC, Associate Professor, Department of Medicine and Division of Cardiology, University of Ottawa, The Ottawa Hospital, Ottawa, Ontario, Canada

#### Correspondence to

Dr Gabriel M Khan;  
DrKhanBooks@gmail.com

#### ABSTRACT

**Objective:** Worldwide, more than 50 million ECGs are carried out annually. Studies show that the standard wrist–ankle placement method may produce poor quality ECGs. Factors including limb movement, tremor, anxiety, cold extremities cause artefacts (distortions of the baseline and ECG waveforms). A problem exists. Physicians face a challenge interpreting poor quality ECGs. Diagnostic delays occur. Patients are subjected to recalls; compassion initiated this study. Torso (modified) leads are favoured as they provide rapid acquisition of better quality ECGs, but cause erroneous electrocardiography. Investigators have not attempted to correct these errors. This study analyses the errors and aims to rectify them.

**Methods:** Our pilot and other studies indicated that torso electrodes placed close to the heart cause an increase in R wave amplitude >3 mm in inferior leads that results in disappearance of inferior infarcts; a decrease >3 mm in I and aVL may cause false lateral infarcts. We hypothesised that finding electrode placements that do not cause alterations in R wave amplitude should correct these problems. Several electrode placements were assessed to obtain ECGs identical with the standard, but with better quality. A total of 1112 patients received standard and new electrode placement recordings. ECG parameters were assessed. A blinded interpreter assisted.

**Results:** Electrodes positioned on the mid-arm and lower abdomen revealed ECGs identical with the standard without artefacts or loss of inferior or appearance of lateral infarcts.

**Conclusions:** This study proved that the new method is faultless and provides better quality ECGs, no recalls and without risk for misinterpretations. Not having to remove leg garments is convenient and allows more rapid acquisition of ECGs. The forearms are freed for, intravenous, radial access, and ECGs needed during procedures. The findings have implications worldwide for patients, clinicians and technicians, and generate the need to study the method in emergency settings.

#### INTRODUCTION

For more than 65 years, the standard 12-lead ECG has been obtained by placing electrodes on the ankles and wrists, along with chest leads<sup>1</sup> and it is well established that limb movement, subtle muscle tremor and interference from muscle potential results in

#### KEY MESSAGES

##### What is already known about this subject?

- It is well established that the Standard 12-lead ECG obtained using the wrist and ankle placement is subject to artefacts (distortions) of the baseline and ECG wave forms caused by limb movement and interference from muscle potential, resulting in poor quality ECGs and thus, the unfortunate recall of patients. Physicians face a major challenge interpreting poor quality ECGs.
- Torso (modified) leads are favoured as they provide rapid acquisition of better quality ECGs.
- Inconvenience of removing hosiery or other leg garments is avoided. However, torso leads cause disappearance of inferior infarcts and the appearance of false lateral infarcts that disqualifies modified leads as a standard recording.

##### What does this study add?

- A solution for torso (modified) leads. This study tested a partial torso lead system that proved the new electrode placement method can replace the wrist–ankle placements for routine electrocardiography.
- Better quality ECGs allow for more precise interpretation, less diagnostic delays, with no recall of patients.
- Compassion for the elderly and infirm who are recalled for repeat ECGs.

##### How might this impact on clinical practice?

- The ECG, a frequently used diagnostic test, is made more reliable. Physicians can provide more precise interpretation for all ECGs. Patients obtain better care.
- The findings generate a reason to initiate studies in emergency settings where lives can be saved worldwide by more rapid acquisition of better quality ECGs and diagnosis that is needed for timely percutaneous coronary intervention or thrombolysis.<sup>2</sup>

artefacts<sup>2, 3</sup> that alter interpretation. Torso (modified) leads have been shown to abolish such artefacts<sup>2, 3</sup> but produce erroneous ECG that disqualifies modified leads as a standard recording. A methodological question remains unresolved. Could modification of the electrode placements used in studies of

modified leads prevent the errors observed? Sound methodology to rectify the problems caused by torso leads required a thorough review and understanding of relevant study findings accompanied by a pilot study.

- ▶ Takuma *et al*<sup>2</sup> in a study of 30 patients, placed electrodes on the shoulder tip and the anterior superior iliac spine to permit easier and more rapid ECG acquisition in emergency settings and obtain higher quality recordings. The electrode placements on the acromion process, and near the iliac crest caused abnormal changes in R wave amplitude. The investigators obtained better quality ECGs and emphasised the advantages of modified leads, but did not recommend them for routine use.<sup>2</sup>
- ▶ Jowett *et al*<sup>3</sup> indicated that ‘The appeal of this easily applied lead system has resulted in extensive use in our hospital, both in emergency and non-emergency situations’.<sup>3</sup> Their study of 100 patients resulted in the disappearance of five of six inferior infarcts that were previously observed on the standard ECG.<sup>3</sup> The electrode sites were as performed by Takuma *et al*.<sup>2</sup> The ECGs of 43% of patients showed movement artefact with wrist–ankle placement compared with no artefacts with torso recordings. R wave amplitude was altered.

Papouchado *et al*<sup>4</sup> assessed the differences between the standard 12-lead ECG and the modified (Mason and Likar<sup>5</sup>) exercise lead system in 29 patients. This showed that when the lower limb electrode was positioned in the nipple line at the level of the umbilicus, and moved vertically to the costal margin or higher, a marked increase in R wave amplitude occurred in inferior leads and a decrease in R wave amplitude was observed in lead I and aVL occurred. Vertical but not horizontal abdominal electrode placements were assessed.<sup>4</sup> Most important, abdominal, horizontal variations were assessed in the present study.

Wilson<sup>6</sup> discussing the standard limb leads placed on the chest stated “the principles of the equilateral triangle<sup>7</sup> should not be applied to any system of chest leads whatsoever”.<sup>6</sup> Wilson emphasised that since all points on the arm have the same potential, it does not matter whether the electrode is placed on the wrist, the forearm or the upper arm; the arm acts mainly as an extension of the wire attached to it.<sup>6</sup>

## METHODS

The acromion process is part of the scapular and thus, part of the trunk. It is not surprising therefore, that this torso electrode position used by Jowett *et al*<sup>3</sup> caused the disappearance of five of six inferior infarcts and appearance of false lateral infarcts and showed significant variation of waveforms in the study by Takuma *et al*.<sup>2</sup> The present study avoided this erroneous electrode placement. Findings from our pilot study, and data from investigators<sup>2–4, 6</sup> initiated the following thinking and deductions to design an

appropriate method to resolve the problems attached to torso (modified) leads.

1. The problem that surrounds torso placements is an increase in R wave amplitude >2 mm in inferior leads that can result in disappearance of inferior myocardial infarction (MI), and a >2 mm decrease in leads I and aVL that may simulate a lateral infarct. Right axis change is caused by alterations in R wave amplitude.
2. Finding electrode placements that do not cause alterations in R wave amplitude should correct the problems.
3. The electrodes for the limb leads must not all be placed on the torso, because this invalidates Einthoven’s<sup>7</sup> theory of the equilateral triangle.<sup>6</sup>
4. Electrodes must not be positioned on the acromion.
5. The upper limb electrode must be placed on the arm at any point<sup>6</sup> but must be below the shoulder joint. This study used the mid-arm; thus, the new electrode placement (NEP) cannot be described as torso leads.

The new electrode placements must prevent the following errors observed with torso leads:


- ▶ An unfavourable change in R wave amplitude in the limb leads that causes:
  - The disappearance of inferior MI;
  - The appearance of QS in leads I and aVL, causing a false lateral infarct;
  - Right axis deviation.

Our pilot study placed electrodes on the mid-arm, various abdominal sites and infraclavicular. Finally, the lower limb electrodes were moved to the lower abdomen 7.6 cm (3 inches) below a horizontal line running across the umbilicus. This position, well below the umbilicus, was maintained for all 1112 patients tested (group A, B, C).

Thought was given to performing ECGs simultaneously, but it was decided that comparing a standard ECG with an NEP ECG performed 1 min apart would not change the interpretation except for arrhythmias. Special software provided duplicate ECGs to be performed at nine laboratory locations, which allowed two electronic ECGs to reach the primary investigator for interpretation. The interpreter was not blinded for the electrodes. ECGs from 1112 patients were assessed. Patients were allocated to group A (786 patients) including 30 with inferior MI, group B (50 patients) and group C (276 patients) including 10 patients with inferior MI.

*Group A:* Enrollment of patients age 33–87 was started following ethics approval, 3 February 2014 in the city of Ottawa. A consent form was signed. Patients for NEP recordings had the ankle electrodes placed 7.6 cm below the umbilicus and 10 cm on either side of the umbilical midline. 651 patients attending Gamma Dynacare laboratories (GDML) for routine ECGs had a standard 12-lead ECG completed; and within 1 min, an NEP ECG was also carried out. A total of 135 patients seen in a cardiology referral centre had duplicate ECGs performed, and were supervised by the primary investigator (figure 1).

**Figure 1** The arm electrodes (A and B) are placed on the mid-arm, on the lateral aspect of the biceps, immediately below the V4 horizontal line. The abdominal electrodes (C and D) placed 7.6 cm (~3 inches) below the umbilical horizontal line, and 5 cm (~2 inches) on either side of the umbilical vertical line. The distance between these two electrodes to be 10 cm (~4 inches).

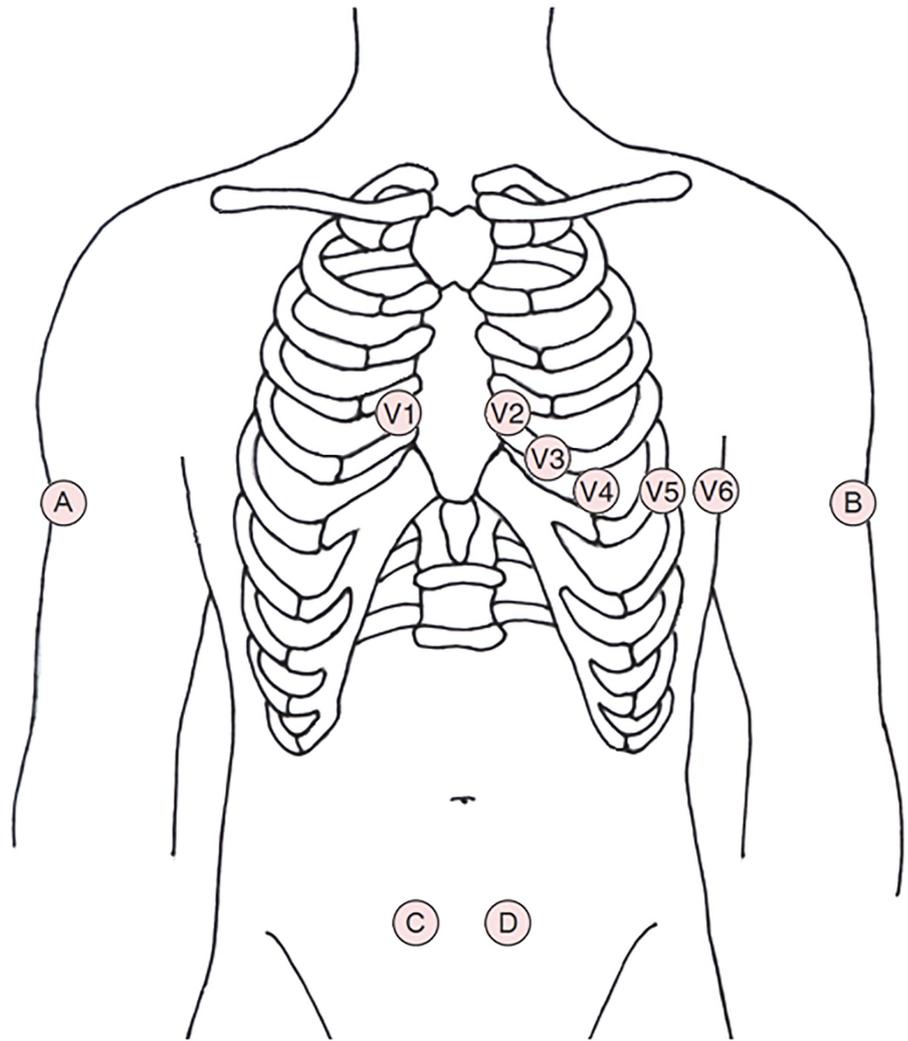

*Duplicate ECGs were assessed for four main problems:*

1. *Inferior MI*—Assessment was made for patients with inferior MIs proven clinically by history of chest pain, serial ECGs, cardiac enzyme elevation, and/or coronary angiograms. One patient was encountered with acute chest pain: the NEP ECG revealed an acute inferior MI; the patient was rushed to coronary intervention; a standard ECG was not completed. For each subject, a standard and NEP ECG was placed side by side on a computer screen and calipers determined if R wave amplitude was increased 2 mm or more in inferior leads. If an old inferior infarct ECG pattern was noted to disappear; or if the R wave amplitude increased 3 mm or more in inferior leads, then the ECG must be repeated under supervision. If the abnormal changes persisted, the study would be halted.
2. *Leads I and aVL*—Assessment was made for a decrease in R wave amplitude of 2 mm or more. A 3 mm or more change may result in a false lateral MI pattern and would halt the study.
3. *Right axis deviation*—The electronic ECGs from >50 consecutive patients in group A were assessed for a greater than 10° right axis change.

4. *Artefacts*—recalls because of artefacts were assessed and poor quality noted.

## RESULTS

Figures 2–5 show that the five major problems observed with torso leads are not observed with recordings carried out using the new placements.


### Group A

1. *Artefacts*—figure 2A shows a standard ECG with artefacts that disappear on the NEP tracing figure 2B (see online supplementary file appendix 3, eFigure 23)

Similar or less prominent artefacts on standard 12-lead ECGs were responsible for 19 (2.94%) recalls of 651 patients but no recalls were required for NEP ECGs. The majority of patient ECGs that revealed major artefacts exhibited clearing of most artefacts when NEP ECGs were completed. Poor quality ECGs were observed with STD ECGs in ~10% of all cases versus ~<2% for NEP ECGs.

2. *Inferior infarcts*—duplicate ECGs from 30 group A patients with proven inferior MI revealed no differences and were reassessed by a blinded interpreter who agreed with our findings. Figure 3 shows the ECG of a patient with a proven old inferior infarct. Diagnostic wide Q

**Figure 2** (A) Standard 12-lead ECG showing artefacts (see online supplementary file appendix 3, eFigure 23). (B). New electrode placement ECG, carried out within 30 s of the standard tracing reveals no artefacts (see online supplementary file appendix 3, eFigure 24).

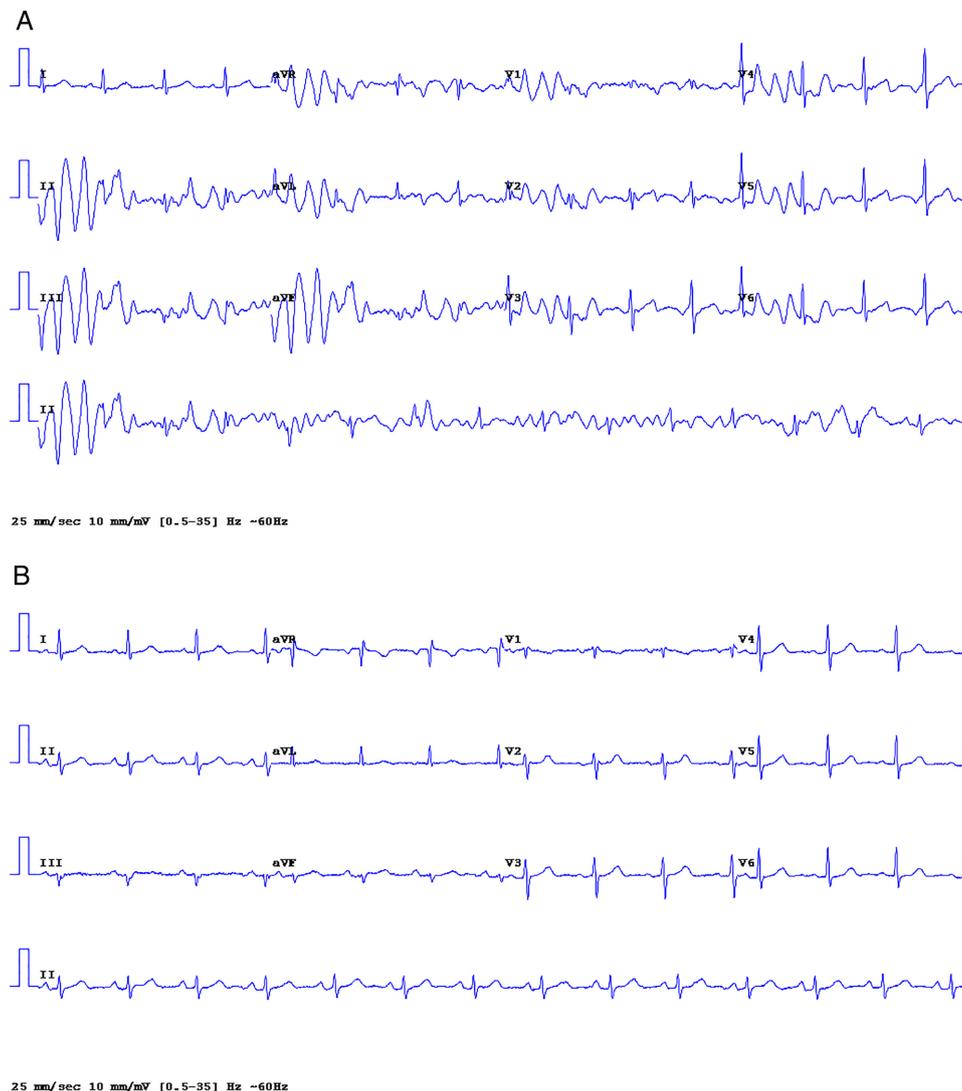

waves are seen in leads 11, 111, aVF: (a) standard ECG, (b) NEP ECG. The NEP ECG of a 57-year-old female patient, with acute chest pain and acute inferior MI, is shown in figure 4. No standard ECG was performed, in order to avoid delays to percutaneous coronary intervention (PCI).


3. *Decrease in R wave amplitude in I and aVL*—No significant decrease in R wave amplitude (equal to or >2 mm) was observed in any ECG from group A (786 patients), and no QS pattern was observed in leads I and aVL. Figure 5A shows a standard ECG tracing with low amplitude (<3 mm) R wave in lead I and aVL. Figure 5B gives the NEP ECG showing no appreciable changes. No QS pattern, in leads I and aVL, false lateral infarct pattern, was observed from the screening of all patient ECGs.

4. *R wave amplitude changes in V leads*: no changes were observed, and no differences in ECGs were noted in patients with right and left bundle branch block, left ventricular hypertrophy or anterior infarction.

5. *Right axis deviation*—fifty group A patients showed a non-significant axis change of zero to 7°

*Group B*—As demonstrated by the pilot study, with abdominal electrode placements as follows: 7.6 cm beneath the umbilicus and horizontally, close to the iliac crest, 20 cm or more from the umbilicus, caused a 1–2 mm change in R wave amplitude versus no change with electrodes 10 cm from the umbilicus. Electrode placements 10 cm from the umbilicus produced recordings that were virtually the same as the standard, and was used successfully in all 786 group A patients. It was felt necessary however, to assess electrodes placed 5 cm from the umbilicus (figure 1). Therefore, on completion of group A, 50 patients (group B) were tested with triplicate ECGs, using placements 10 cm and 5 cm from the umbilicus compared with standard recordings. For each sets of recordings, the technician was supervised by the investigator.

*Results Group B*: the abdominal placement 10 cm from the umbilicus showed ECGs virtually identical with the standard in all patients. The R wave amplitude was equal to the standard in 44 of 50 patients. In six patients minimal, 0.5–1 mm non-significant R wave changes were observed. Electrode placement 5 cm from the umbilicus

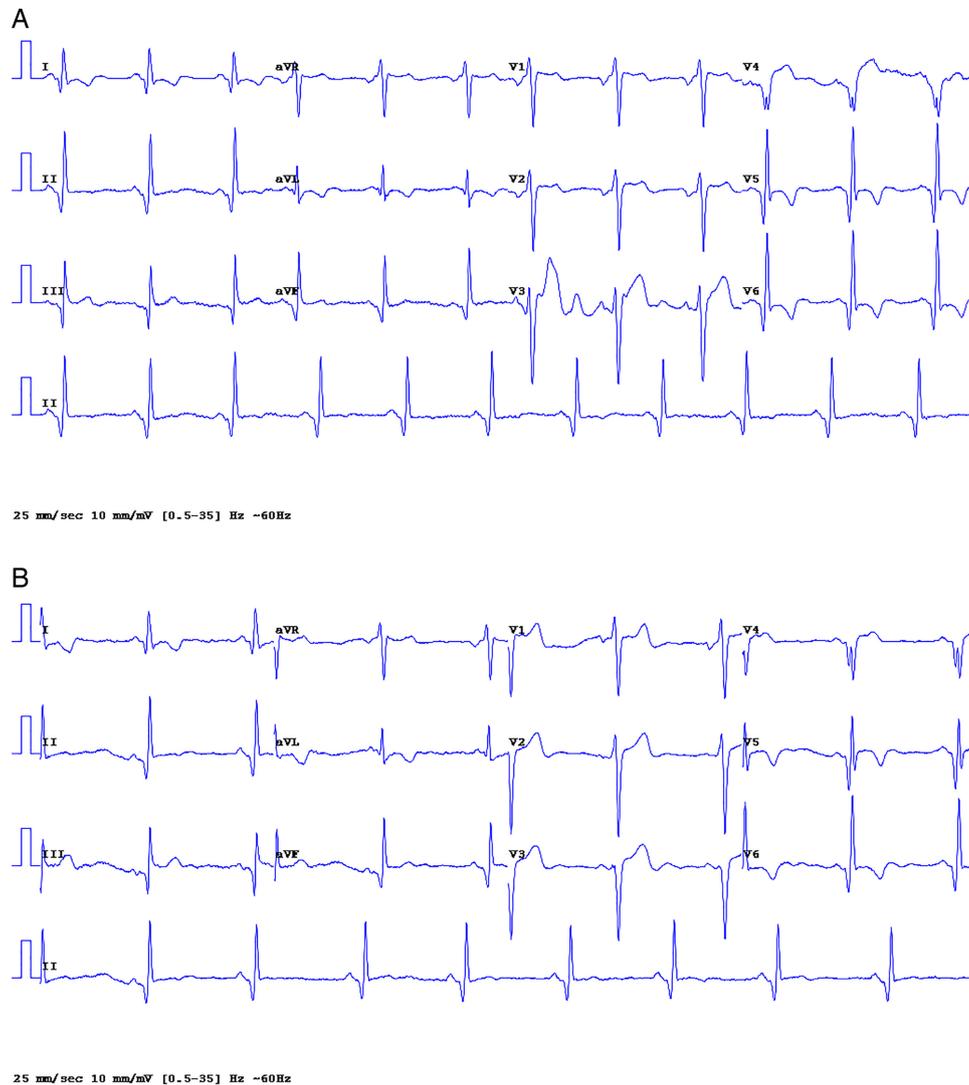

**Figure 3** (A) Standard lead ECG: deep wide Q waves in leads 11, 111, aVF: old inferior, and anterolateral myocardial infarction, artefacts in V3, V4 (see online supplementary file, appendix 1, eFigures 1–59, eAppendix 3, figures 1–20). (B) New electrode placement ECG, old inferolateral infarct; the tracing is similar to the standard recording, but with clearing of artefacts in V3 and V4. The R wave amplitude in all 12 leads is the same in the standard and NEP ECG.

(figure 1) revealed identical wave forms and R wave amplitude as the standard in all 50. Although R wave amplitude changes at 10 cm from the umbilicus were considered not significant, and was used successfully in all group A patients, group C was started with abdominal placement 5 cm from the umbilicus (figure 1). Group B patients were, in essence, a second pilot study and were not counted in the statistical analysis.


*Group C*—276 patients attending GDML had duplicate ECGs performed with electrode placements as shown in figure 1.

*Results* The principal investigator, and a blinded interpreter compared NEP and standard ECGs and found no differences in wave forms, or parameters. Table 1 gives ECG parameters in .30 consecutive patients. No changes in R wave amplitude were noted, and no false lateral infarcts appeared. The ECGs of 10 patients with proven old inferior MIs revealed the same waveforms, with no

disappearance of inferior infarcts (see online supplementary file appendix 3). P waves, ST segment and T waves were identical with that of the standard ECG. No ECG differences were noted in patients with left and right bundle branch block, atrial fibrillation or premature beats. Arrhythmia diagnosis was unchanged by the new method. No R wave amplitude changes occurred in leads 1 and aVL, inferior or chest leads. Forty-six patients had low amplitude R waves <5 mm in leads 1 and aVL, and 16 with the lowest R wave amplitudes <4 mm were chosen to test the new placement method for flaws. No occurrence of a QS pattern in leads 1 and aVL or false lateral infarct was noted in the 16 test patients ECGs (see test cases in the online supplementary file, appendix 2, eFigures 1 to 32 and appendix 3 eFigure 21, 22).


A blinded interpreter assessed:

- ▶ All 276 group C patients.

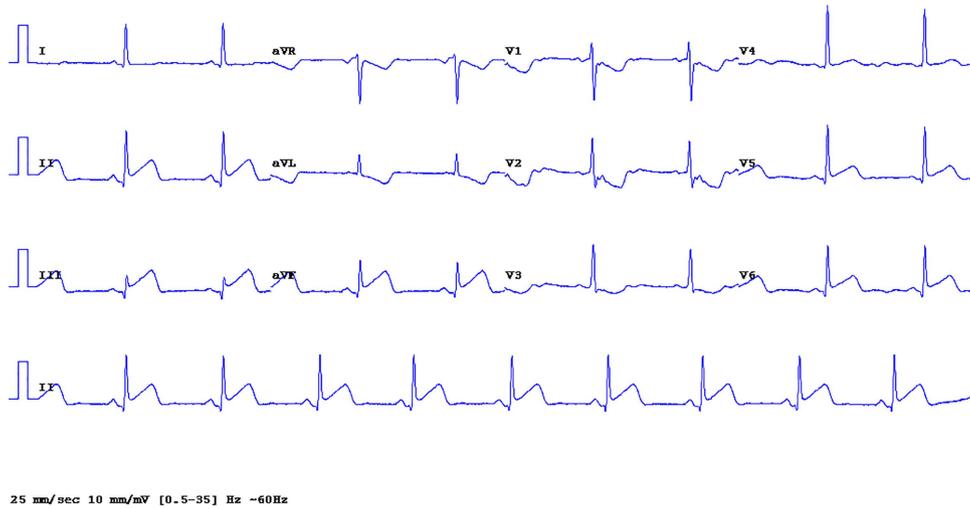

**Figure 4** New electrode placement ECG from a female, age 57 with acute chest pain, shows abnormal ST elevation in leads I, II, III, aVF, V5, V6: typical findings of acute inferior MI with lateral involvement, ST segment elevation MI. Reciprocal depression V1-V2, aVL is not diagnostic but helps confirm the diagnosis of acute MI.

**Figure 5** (A) Standard 12-lead ECG: low R wave amplitude in lead I, aVL; a test case chosen to flaw the new electrode placement (see online supplementary file, appendix 2, eFigures 1 to 32 and appendix 3 eFigure 21, 22). (B) New electrode placement ECG; similar to that of the standard ECG. The R wave amplitude in all 12 leads is virtually the same. The very small R wave in aVL did not change to a QS pattern (see online supplementary file eAppendix 2, eFigures 1 to 32).

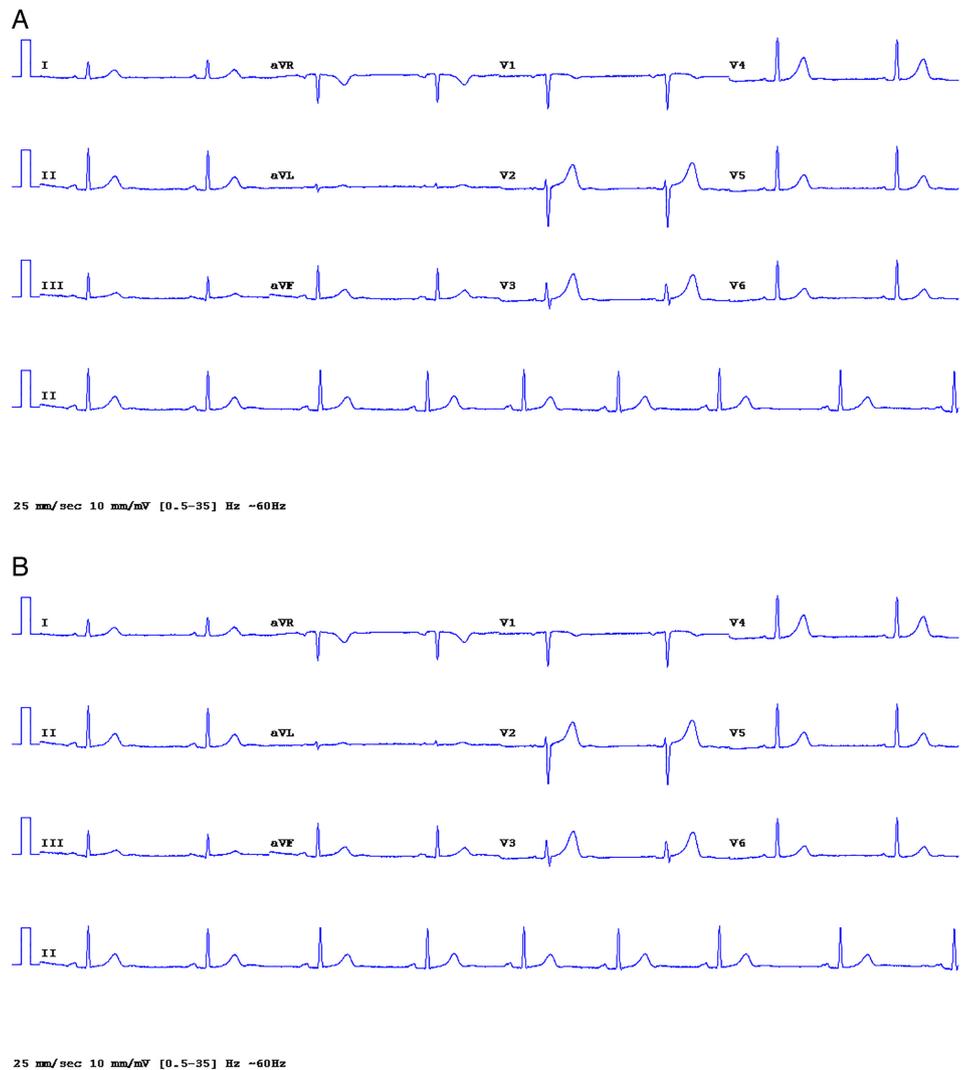

**Table 1** R wave amplitude in lead 1 and ECG parameters in 30 patients**R wave amplitude and ECG parameters in 30 patients from group C (276 patients)**

Subject Parameters/ECG type	1		2		3		4		5	
	NEP	STD								
PR interval (ms)	177	178	152	155	172	173	173	172	148	150
QRS duration	98	95	98	100	110	112	105	101	100	92
Axis	-17	-19	-3	-7	-31	-27	-11	-8	-38	-33
R wave amplitude lead 1 (mm)	14	14	10	11	6	5	13	15	6	6
QTc	404	402	438	424	443	443	400	406	435	427
Subject	6		7		8		9		10	
PR interval (ms)	200	200	130	128	158	148	170	168	173	170
QRS duration	132	133	82	87	67	67	110	115	129	132
Axis	-83	-84	55	48	83	82	65	67	89	91
R wave amplitude lead 1 (mm)	3	3.5	7	7	5	5	5	5	4	4
QTc	432	430	392	398	382	385	444	442	452	448
Subject	11		12		13		14		15	
PR interval (ms)	158	160	225	227	203	197	150	143	173	177
QRS duration	84	83	175	173	120	135	93	97	180	190
Axis	4	14	-65	-63	-41	-40	6	4	82	81
R wave amplitude lead 1 (mm)	10	10	3	3	11	11	7	7	3	3
QTc	402	403	496	482	425	452	396	399	499	505
Subject	16		17		18		19		20	
PR interval (ms)	150	153	168	155	130	133	155	153	193	187
QRS duration	97	97	105	118	183	197	85	84	114	113
Axis	36	39	-13	-15	13	6	68	64	62	55
R wave amplitude lead 1 (mm)	7	6	7	7	13	13	4	4.5	6.5	7
QTc	400	409	419	424	499	493	412	411	420	425
Subject	21		22		23		24		25	
PR interval (ms)	187	188	143	140	158	162	147	137	187	210
QRS duration	130	129	97	100	87	85	130	140	93	97
Axis	-24	-25	25	22	39	40	-54	-49	-24	-21
R wave amplitude lead 1 (mm)	10	9	9	9	5	5	7	7	7	7
QTc	418	409	406	435	433	436	448	445	434	440
Subject	26		27		28		29		30	
PR interval (ms)	183	180	173	173	157	147	150	147	223	220
QRS duration	87	90	124	114	76	86	93	86	90	90
Axis	68	68	40	44	23	18	51	50	37	34
R wave amplitude lead 1 (mm)	4	4	3	3	8	8	7	7.5	6.5	6.5
QTc	430	411	434	434	385	392	431	432	400	408

- ▶ All ECGs that are in the online supplementary files.
  - 10 duplicate ECGs with inferior MIs (see online supplementary file appendix 3).
- ▶ The QRS axis in 50 consecutive group C patients (see online supplementary file table S2).
  - 16 duplicate ECGs, test cases showing very low R waves in leads 1 and aVL for QS pattern or false lateral MI (see online supplementary file eAppendix 2, eFigures 1 to 32).
  - 30 ECGs showing old inferior MIs (see online supplementary file eAppendix 1).

Figure 6A shows the standard 12-lead ECG from a patient with Eisenmenger syndrome, marked right ventricular hypertrophy: R >1 mV in V1, R/S >1 in V1, R/S <1 in V4, V5, V6, right atrial hypertrophy and right axis deviation 153. The NEP ECG shown in figure 6B reveals virtually the same waveforms and QRS axis of 155.

### STATISTICAL ANALYSIS

Intraclass correlation (ICC) estimates were calculated in order to quantify the agreement between the NEP ECG

and the standard (STD) ECG with respect to R wave amplitude in Lead 1 and the different ECG parameters measured in n=30 patients. The NEP and STD method have a high level of agreement in measurement of the R wave amplitude in Lead 1, ICC correlation 0.984. **Table 1** shows R wave amplitudes in Lead 1, and ECG parameters for the NEP and standard ECGs in 30 patients. See **table 2** for ICC estimates and Spearman correlations. Spearman and Pearson correlations as well as the ICC were also calculated in order to quantify the agreement in the axis deviation measured in 50 consecutive patients (See **table 2A**, **table 1** in the supplement gives actual measurements of axis deviation).

When comparing diagnosis with STD and NEP, sensitivities and specificities were also calculated along with exact binomial 95% CIs (**table 3**). The NEP and standard ECG recordings from 40 patients with inferior infarction and 149 patients with normal ECGs show similar sensitivities and specificities (**table 3**).

## DISCUSSION

This large study proved conclusively that the NEP method is sound and suitable for routine ECG. Five small studies compared standard limb lead placements and torso placed electrodes. Diamond *et al*<sup>8</sup> studied only 11 patients using all torso leads. Variations in QRS amplitude were observed but reported as not significant. The investigators concluded that the torso lead system appears to produce an ECG recording interchangeable with the standard limb placement. Time has shown that this conclusion is incorrect. Importantly, there were no patients with inferior infarction in the study. The standard ECG in one patient 'with a marked tremor was interpretable only with great difficulty; the torso system gave a far superior ECG'.<sup>8</sup>

A 1989 study by Sevilla *et al*<sup>9</sup> assessed 44 patients, using exercise torso-lead placements to acquire 12-lead ECGs. The wrist electrodes were placed under the clavicle in the mid-clavicular lines. The ankle electrodes

were positioned in the anterior axillary lines above the costal margins. Loss of ECG inferior pattern occurred in 11 of 16; a false lateral infarct was seen in 19% of patients and a significant change in frontal axis was noted.

Pahlm *et al*<sup>10</sup> reported on torso placements in 30 patients. Simultaneous recordings revealed significant frontal plane axis shift and diminished Q-wave durations in lead aVF compared with those of the standard.

The 1995 emergency room study of 30 patients with only one inferior infarct by Takuma *et al*<sup>2</sup> showed minor QRS amplitude changes, and was not recommended for routine use. However, the quality of recordings were definitely improved.<sup>2</sup> The wrist electrodes were placed on the shoulder tip, the acromion process.

The 2005, ECG study of 100 patients by Jowett *et al*<sup>3</sup> with torso electrode similar to the Takuma placement resulted in loss of five of six inferior infarcts, one false lateral infarct, and important amplitude changes that created eight possible myocardial infarcts to appear. A marked change in axis was observed.

Importantly, these investigators<sup>2 3 9 10</sup> did not attempt to correct the problems caused by torso leads. The present study analysed the problems and rectified them.

The 2007 AHA/ACC guidelines on ECG Part 1: recommends studies to look at modified leads.<sup>11</sup> Previous studies<sup>3 4 9</sup> and our pilot study indicate that modified leads (all torso) cause serious errors.

The present study shows that a lower torso and upper arm placement is reliable. No inferior infarcts were lost and no false lateral infarcts were produced by the NEP. If the abdominal electrode placement is as shown in **figure 1**, the NEP does not cause a change in R wave amplitude. Thus, disappearance of inferior MI and causation of a lateral infarct cannot occur. Each ECG should be labelled 'new electrode placement' to distinguish the system from modified leads that provide erroneous ECG recordings. However, a 5 cm lateral displacement of the abdominal electrodes should not alter interpretation as shown for all 786 group A patients including 30 patients with inferior infarcts in whom the abdominal electrode placement was ~5 cm lateral to that used for group C (see **figure 1**).

## Limitations and further studies

This study was not performed in emergency rooms, but the NEP ECG captured the ST-elevation, reciprocal depression and Q waves of acute inferior MI remarkably well (**figure 4**). Owing to the rapidity of application, and better quality ECGs, the NEP should prove useful for diagnosis in patients with probable acute coronary syndrome and in emergency settings. The NEP needs to be studied in this population of patients.

## CONCLUSION

This large study showed that the new method produces ECGs that are identical to the standard and provides

**Table 2** Intraclass correlation (ICC) and Spearman correlation in 30 patients showing R wave amplitude and ECG parameters and axis deviation in 50 patients

	ICC intraclass correlation	Pearson	Spearman
ECG parameters (n=30 patients)			
PR interval	0.965		0.967
QRS duration	0.977		0.953
QRS axis	0.996		0.996
R wave amplitude lead 1	0.984		0.982
QTc	0.943		0.899
Axis deviation (n=50 consecutive patients)	0.992	0.993	0.990



**Figure 6** (A) Standard 12-lead ECG from a patient with Eisenmenger syndrome. Marked right ventricular hypertrophy:  $R > 1$  mV in V1,  $R/S > 1$  in V1,  $R/S < 1$  in V4, V5, V6, right atrial hypertrophy, and right axis deviation 153. (B) New electrode placement ECG from the patient shown in figure 6A, QRS axis 155; findings are identical to the standard.

better quality ECGs, with no recalls and without risk for misinterpretations.

- ▶ The NEP revealed waveforms that are identical with that of the standard 12-lead ECG, resulting in the same interpretation;
- ▶ The new method does not cause disappearance of MI pattern or change R wave amplitude or QRS axis;

- ▶ The NEP was shown suitable and as reliable as the standard for observation of P waves and arrhythmia diagnosis;
- ▶ The NEP can replace standard electrode placements at the wrist and ankles for routine 12-lead electrocardiography.

The NEP is an alternative method to the wrist–ankle placements and has several advantages:

**Table 3** Sensitivity and specificity of the NEP method for diagnosis of old inferior MI and normal ECG

Diagnosis					Sensitivity (95% CI)	Specificity (95% CI)
Old inferior MI (n=1112)					1.0 (0.997 to 1.0)	1.0 (0.912 to 1.0)
Normal ECG (group C) (n=276)					1.0 (0.971 to 1.0)	1.0 (0.976 to 1.0)
Diagnosis	TP	FP	TN	FN	Sensitivity (95% CI)	Specificity (95% CI)
Old inferior MI (n=1112)	40	0	1072	0	1.0 (0.997 to 1.0)	1.0 (0.912 to 1.0)
Normal ECG (group C) (n=276)	149	0	127	0	1.0 (0.971 to 1.0)	1.0 (0.976 to 1.0)

FN, False negative; FP, False positive; TN, True negative; TP, True positive.

- ▶ Recordings with better quality, no recalls;
- ▶ Not having to remove hosiery or other leg garments is convenient and allows more rapid acquisition of ECGs. Time saved is important for therapy, and for hospital and non-hospital laboratories.

The findings have implications worldwide for patients, clinicians, technicians and medical laboratories, and generate the need to study the NEP in patients with probable acute coronary syndrome where timely PCI or thrombolysis can save lives. In addition, the forearms are freed for intravenous or radial artery access and further ECGs that are needed during procedures can be obtained.

**Acknowledgements** The study could not have enrolled this large number of patients, and completion of Group C could not be accomplished without the help of Lanor Hitch, MLT Manager of Cardiology Services Gamma Dynacare Medical Laboratories. Lanor and her supervisors accomplished the training of technicians for the study. Professor George A Wells, School of Epidemiology, Public Health and Preventative Medicine, University of Ottawa. Director of the Cardiovascular Research Methods Centre, University of Ottawa Heart Institute, and Jordan Bernick MSc provided the statistical analysis. Validation of the study results were completed by an independent interpreter, unknown to the study investigator: Dr Timothy Burns MD, FRCP C, diligently assessed all ECGs from group C patients, along with 40 ECGs showing inferior infarcts, and all shown in the supplement without obtaining remuneration or gifts. Dr Michele A Turek, Associate Professor, Department of Medicine and Division of Cardiology, University of Ottawa, provided editorial and content support. Stephen Khan assisted with editing, revisions, reformatting of the manuscript, tables and indexing the supplementary material. The author would like to thank Dr Max della Zazzera Medical Director, Bank Medical Centre, for allowing him to use a technician to complete triplicate ECG recordings on patients for the study. Yasmin Khan was responsible for the art work (figure 1). Carrie Barlow assisted with editing of the initial manuscript.

**Contributors** The Study, design, acquisition, interpretation of data, and writing the manuscript were carried out by GMK who had full access to all of the data in the study and takes responsibility for the integrity of the data and the accuracy of the data analysis.

**Competing interests** None declared.

**Patient consent** Obtained.

**Ethics approval** Gamma-Dynacare Medical Laboratories (GDML).

**Provenance and peer review** Not commissioned; externally peer reviewed.

**Data sharing statement** Article Title: A New Electrode Placement Method for obtaining 12-lead ECGs Author: GMK M.D. Email for further information: DrKhanBooks@gmail.com Unpublished Data: 1. ECG artefacts are a major

problem in out of hospital laboratories 2. Modified lead placement for Parkinson's disease, tremors, amputees, and patients who must remain seated. 3. Cold extremities and artefacts 4. How big a problem do we face.

**Open Access** This is an Open Access article distributed in accordance with the Creative Commons Attribution Non Commercial (CC BY-NC 4.0) license, which permits others to distribute, remix, adapt, build upon this work non-commercially, and license their derivative works on different terms, provided the original work is properly cited and the use is non-commercial. See: <http://creativecommons.org/licenses/by-nc/4.0/>

## REFERENCES

1. Wilson FN, Johnson FD, Rosenbaum FF, *et al.* The precordial electrocardiogram. *Am Heart J* 1944;27:19–85.
2. Takuma K, Hori S, Sasaki J, *et al.* An alternative limb lead system for electrocardiographs in emergency patients. *Am J Emerg Med* 1995;13:514–17.
3. Jowett NI, Turner AM, Cole A, *et al.* Modified electrode placement must be recorded when performing 12-lead electrocardiograms. *Postgrad Med J* 2005;81:122–5.
4. Papouchado M, Walker PR, James MA, *et al.* Fundamental differences between the standard 12-lead electrocardiograph and the modified (Mason-Likar) exercise lead system. *Eur Heart J* 1987;8:725–33.
5. Mason RE, Likar I. A new system of multiple-lead exercise electrocardiography. *Am Heart J* 1966;71:196–205.
6. Wilson FN. The distribution of the potential differences produced by the heart beat within the body and its surface. *Am Heart J* 1929–1930;5:599–616.
7. Einthoven W, Fahr G, de Waart A. On the direction and manifest size of the variations of potential in the human heart and on the influence of the position of the heart on the form of the electrocardiogram. *Pflugers Arch* 1913;150:275–315, [Translation by: Hoff HE, Sekelj P. *Am Heart J* 1950;40(2):163–211].
8. Diamond D, Griffith DH, Greenberg ML, *et al.* Torso mounted electrocardiographic electrodes for routine clinical electrocardiography. *J Electrocardiol* 1979;12:403–4069.
9. Sevilla DC, Dohrmann ML, Somelofski CA, *et al.* Invalidation of the resting electrocardiogram obtained via exercise electrode sites as a standard 12-lead recording. *Am J Cardiol* 1989;63:35–9.
10. Pahlm O, Haisty WK Jr, Edenbrandt L, *et al.* Evaluation of changes in standard electrocardiographic QRS waveforms recorded from activity-compatible proximal limb lead positions. *Am J Cardiol* 1992;69:253–7.
11. Kligfield P, Gettes LS, Bailey JJ, American Heart Association Electrocardiography and Arrhythmias Committee, Council on Clinical Cardiology; American College of Cardiology Foundation; Heart Rhythm Society, Josephson M, Mason JW, Okin P, Surawicz B, Wellens H. Recommendations for the standardization and interpretation of the ECG: part I: The electrocardiogram and its technology: a scientific statement from the American Heart Association Electrocardiography and Arrhythmias Committee, Council on Clinical Cardiology; the American College of Cardiology Foundation; and the Heart Rhythm Society: endorsed by the International Society for Computerized Electrocardiology *Circulation* 2007;115:1306–24.

## SUPPLEMENTAL MATERIAL

Khan, Gabriel M.

A New Electrode Placement Method for 12-lead Electrocardiography.

The supplement gives standard (STD) and new electrode placement (NEP) ECGs that were verified by a blinded interpreter. They focus on the four errors known to occur with torso (modified) leads).

Marked changes in R wave amplitude causing:

- disappearance of inferior infarcts;
- appearance of QS in leads 1, aVL – probable false lateral infarcts;
- right axis deviation.

eAppendix 1. 30 ECGs from patients in group A, all with proven inferior MI. The abdominal electrode was placed 3 inches below the umbilical level and 8 inches apart. Standard and new electrode placement ECGs show identical wave forms, and similar interpretations.

eAppendix 2. 16 ECGs from patients in group A: these ECGs were chosen from 786 group A patient's ECGs because of the presence of very low R waves in leads 1 and aVL; that may reveal a QS on the NEP ECG simulating a false lateral MI; test cases to fault the new electrode placement method.

eAppendix 3. 13 ECGs from patients in group C: 10 with inferior MI, 1 with low R waves lead 1, AVL ; 1 patient with severe right ventricular hypertrophy and marked right axis deviation; 1 patient with standard ECG showing artifacts that clear on NEP ECG, a better quality ECG.

All group C patients 1 to 264 and all ECGs provided in appendix 1 to 3 have been validated by a blinded interpreter unknown to the author.

**Table 1.** QRS Axis deviation in 50 consecutive patients

**eAppendix 1.** 30 subjects with proven inferior MI: test cases from group A; abdominal placement 8 inches apart: Standard and new electrode placement ECGs are virtually identical. All ECGs verified by blinded- interpreter.

**eFigure 1.** Patient 01 old Inferior MI STD ECG

**eFigure 2.** Patient 01 old Inferior MI NEP ECG

**eFigure 3.** Patient 02 old inferior MI STD ECG

**eFigure 4.** Patient 02 old inferior MI NEP ECG

**eFigure 5.** Patient 03 old inferior MI STD ECG

**eFigure 6.** Patient 03 old inferior MI NEP ECG

**eFigure 7.** Patient 04 old inferior MI STD ECG  
**eFigure 8.** Patient 04 old inferior MI NEP ECG  
**eFigure 9.** Patient 05 old inferior MI STD ECG  
**eFigure 10.** Patient 05 old inferior MI NEP ECG  
**eFigure 11.** Patient 06 old inferior MI STD ECG  
**eFigure 12.** Patient 06 old inferior MI NEP ECG  
**eFigure 13a.** Patient 07 old inferior MI STD ECG  
**eFigure 13b.** Patient 07 old inferior MI NEP ECG  
**eFigure 14.** Patient 08 old inferior MI STD ECG  
**eFigure 15.** Patient 08 old inferior MI NEP ECG  
**eFigure 16.** Patient 09 old inferior MI STD ECG  
**eFigure 17.** Patient 09 old inferior MI NEP ECG  
**eFigure 18.** Patient 10 old inferior MI STD ECG  
**eFigure 19.** Patient 10 old inferior MI NEP ECG  
**eFigure 20.** Patient 11 old inferior MI STD ECG  
**eFigure 21.** Patient 11 old Inferior MI NEP ECG  
**eFigure 22.** Patient 12 old inferior MI STD ECG  
**eFigure 23.** Patient 12 old inferior MI NEP ECG  
**eFigure 24.** Patient 13 old inferior MI STD ECG  
**eFigure 25.** Patient 13 old inferior MI NEP ECG  
**eFigure 26.** Patient 14 old inferior MI STD ECG  
**eFigure 27.** Patient 14 old inferior MI NEP ECG  
**eFigure 28.** Patient 15 old inferior MI STD ECG  
**eFigure 29.** Patient 15 old inferior MI NEP ECG  
**eFigure 30.** Patient 16 old inferior MI STD ECG  
**eFigure 31.** Patient 16 old inferior MI NEP ECG  
**eFigure 32.** Patient 17 old inferior MI NEP ECG  
**eFigure 33.** Patient 17 old inferior MI STD ECG  
**eFigure 34.** Patient 18 old inferior MI STD ECG  
**eFigure 35.** Patient 18 old inferior MI NEP ECG  
**eFigure 36.** Patient 19 old inferior MI STD ECG  
**eFigure 37.** Patient 19 old inferior MI NEP ECG  
**eFigure 38.** Patient 20 old inferior MI STD ECG  
**eFigure 39.** Patient 20 old inferior MI NEP ECG  
**eFigure 40.** Patient 21 old inferior MI STD ECG  
**eFigure 41.** Patient 21 old inferior MI NEP ECG  
**eFigure 42.** Patient 22 old inferior MI STD ECG  
**eFigure 43.** Patient 22 old inferior MI NEP ECG  
**eFigure 44.** Patient 23 old inferior MI STD ECG  
**eFigure 45.** Patient 23 old inferior MI NEP ECG  
**eFigure 46.** Patient 24 old inferior MI STD ECG  
**eFigure 47.** Patient 24 old inferior MI NEP ECG  
**eFigure 48.** Patient 25 inferior MI NEP ECG similar to STD  
**eFigure 49.** Patient 25 proven inferior MI STD ECG  
**eFigure 50.** Patient 26 inferior MI NEP ECG similar to STD  
**eFigure 51.** Patient 26 proven inferior MI STD ECG

- eFigure 52.** Patient 27 old inferior MI STD ECG
- eFigure 53.** Patient 27 old inferior MI NEP ECG
- eFigure 54.** Patient 28 proven inferior MI STD ECG
- eFigure 55.** Patient 28 inferior MI NEP ECG similar to STD
- eFigure 56.** Patient 29 old inferior MI STD ECG
- eFigure 57.** Patient 29 old inferior MI NEP ECG
- eFigure 58.** Patient 30 old inferior MI STD ECG artifacts
- eFigure 59.** Patient 30 old inferior MI NEP ECG artifacts cleared

**eAppendix 2.** 16 test cases from group a: low R wave amplitude in leads 1 and aVL in standard 12 – lead ECG; NEP ECG shows no decrease in R amplitude or conversion to QS pattern.

- eFigure 1.** Patient 01 lead 1 aVL STD ECG R waves similar to NEP
- eFigure 2.** Patient 01 lead 1 aVL NEP ECG R wave amplitude is similar to STD ECG
- eFigure 3.** Patient 02 lead 1 aVL STD ECG similar R wave amplitude to NEP
- eFigure 4.** Patient 02 LEAD 1 AVL NEP ECG, R wave amplitude in leads 1, aVL and other leads is similar to the STD ECG
- eFigure 5.** Patient 03 lead 1 aVL STD ECG very low R wave amplitude leads 1 aVL
- eFigure 6.** Patient 03 lead 1aVL NEP ECG no change in R wave amplitude compared with STD
- eFigure 7.** Patient 04 lead 1 aVL STD ECG
- eFigure 8.** Patient 04 lead 1 aVL NEP ECG
- eFigure 9.** Patient 05 lead 1 VL STD ECG
- eFigure 10.** Patient 05 lead 1 aVL NEP ECG R wave amplitude in all leads are similar to STD
- eFigure 11.** Patient 06 lead 1 aVL STD ECG
- eFigure 12.** Patient 06 lead 1 aVL NEP ECG R wave amplitude in all 12 leads similar to STD
- eFigure 13.** Patient 07 lead 1 aVL NEP ECG all R waves similar to STD ECG
- eFigure 14.** Patient 07 lead 1 aVL STD ECG
- eFigure 15.** Patient 08 lead 1 aVL STD ECG
- eFigure 16.** Patient 08 lead 1 aVL NEP ECG similar to STD
- eFigure 17.** Patient 09 lead 1 aVL STD ECG
- eFigure 18.** Patient 09 lead 1 aVL NEP ECG similar
- eFigure 19.** Patient 10 lead 1 aVL STD ECG
- eFigure 20.** Patient 10 lead 1 aVL NEP ECG similar R waves
- eFigure 21.** Patient 11 LEAD 1 AVL STD ECG
- eFigure 22.** Patient 11 lead 1 aVL NEP ECG all R waves amplitude similar to STD
- eFigure 23.** Patient 12 lead 1 aVL STD ECG
- eFigure 24.** Patient 12 lead 1 aVL NEP ECG
- eFigure 25.** Patient 13 lead 1 aVL STD ECG
- eFigure 26.** Patient 13 lead 1 aVL NEP ECG R wave amplitude in inferior and other leads similar to STD
- eFigure 27.** Patient 14 lead 1 aVL STD ECG
- eFigure 28.** Patient 14 lead 1 aVL NEP ECG similar
- eFigure 29.** Patient 15 lead 1 aVL STD ECG

**eFigure 30.** Patient 15 Patient 15 LEAD 1 AVL NEP ECG virtually identical R wave amplitude all leads.

**eFigure 31.** Patient 16 R lead 1 aVL STD ECG

**eFigure 32.** Patient 16 R wave amplitude NEP ECG no increase in R wave amplitude in 11 111 aVF and no decrease in lead 1 Avl.

**eAppendix 3.** 10 proven old Inferior MI, and 3 test cases of interest from group C: standard and new electrode placement ECGs.

**eFigure 1.** 2014/05/30 Standard ECG, deep, wide inferior Q waves; old inferior infarct, QRS axis 19.

**eFigure 2.** NEP ECG shows virtually the same findings, deep, wide Q waves: old inferior MI; QRS axis 25.

**eFigure 3.** 2014/ 05/09; standard ECG old inferior MI, QRS axis 74. VPBs.

**eFigure 4.** new electrode placement ECG, similar Q waves in inferior leads, axis 69.

**eFigure 5.** standard ECG, 2014/05/23; Q waves 11,111, aVF: old inferior infarct, VPBs.

**eFigure 6.** NEP ECG reveals virtually the same, Q waves 11, 111, aVF, old inferior MI, VPBs.

**eFigure 7.** 2014 /05 /16 Standard 12 – lead ECG, Q waves 11, 111, aVF, old inferior MI

**eFigure 8.** New electrode placement ECG, similar wave form and inferior Q waves

**eFigure 9.** 14/05/23 Standard ECG; old inferior MI, QRS axis -40. RBBB

**eFigure 10.** NEP ECG similar wave form, inferior Q waves, axis -35, RBBB. Virtually the same

**eFigure 11.** 14/05/27; standard ECG; Q waves 11, 111 aVF; old inferior MI, axis -30

**eFigure 12.** NEP ECG similar inferior Q waves, axis -30; R wave amplitude in lead 1 aVL unchanged

**eFigure 13.** Standard ECG, old inferior MI, axis -57 virtually the same wave form as NEP ECG shown below

**eFigure 14.** New electrode placement ECG; similar inferior Q waves; old inferior MI, axis -55

**eFigure 15.** Standard ECG, old inferior MI; axis -47, note artifacts

**eFigure 16.** NEP ECG, old inferior MI similar to standard recording, artifacts abolished, axis - 42

**eFigure 17.** Standard ECG, inferior Q waves: old inferior MI

**eFigure 18.** NEP ECG, old inferior MI, virtually similar wave forms, inferior Q waves

**eFigure 19.** Standard ECG old inferior MI, right coronary occlusion, quadruple CABG

**eFigure 20.** NEP ECG, old inferior MI, similar

**eFigure 21.** Standard 12 – lead ECG; low amplitude R waves in leads 1 and aVL that can result in a QS pattern as observed with torso ( modified ) - leads

**eFigure 22.** New electrode placement ECG showing no change in R wave amplitude in leads 1, aVL or inferior leads. Virtually the same wave forms.

**eFigure 23.** A standard method ECG shows artifacts; poor quality ECG

**eFigure 24.** NEP ECG done within 1 minute of STD shows clearing of artifacts a better quality ECG.

This supplementary material has been provided by the author to give readers additional information about the work.

Axis Deviation (50 consecutive subjects)					
Subject	NEP	STD	Subject	NEP	STD
1	-36	-41	26	-18	-25
2	-2	-2	27	42	35
3	65	64	28	10	7
4	44	44	29	58	54
5	27	31	30	64	57
6	21	22	31	-19	-18
7	-21	-24	32	16	14
8	-25	-24	33	18	10
9	-13	-15	34	82	85
10	83	84	35	-2	-3
11	-17	-10	36	2	-5
12	52	45	37	90	90
13	52	47	38	63	62
14	55	52	39	52	41
15	9	16	40	-1	-4
16	10	0	41	-82	-81
17	-14	-9	42	74	78
18	77	74	43	30	34
19	47	40	44	80	75
20	33	24	45	41	40
21	61	59	46	79	75
22	23	13	47	155	153
23	83	80	48	45	31
24	9	0	49	52	44
25	78	77	50	56	60

**Table 1.** QRS Axis deviation in 50 consecutive patients


## eAppendix 1



eFigure 1: Patient 01; deep wide inferior Q waves; old Inferior MI STD ECG



eFigure 2: Patient 01 ~ 20 year old Inferior MI; the NEP ECG is identical.

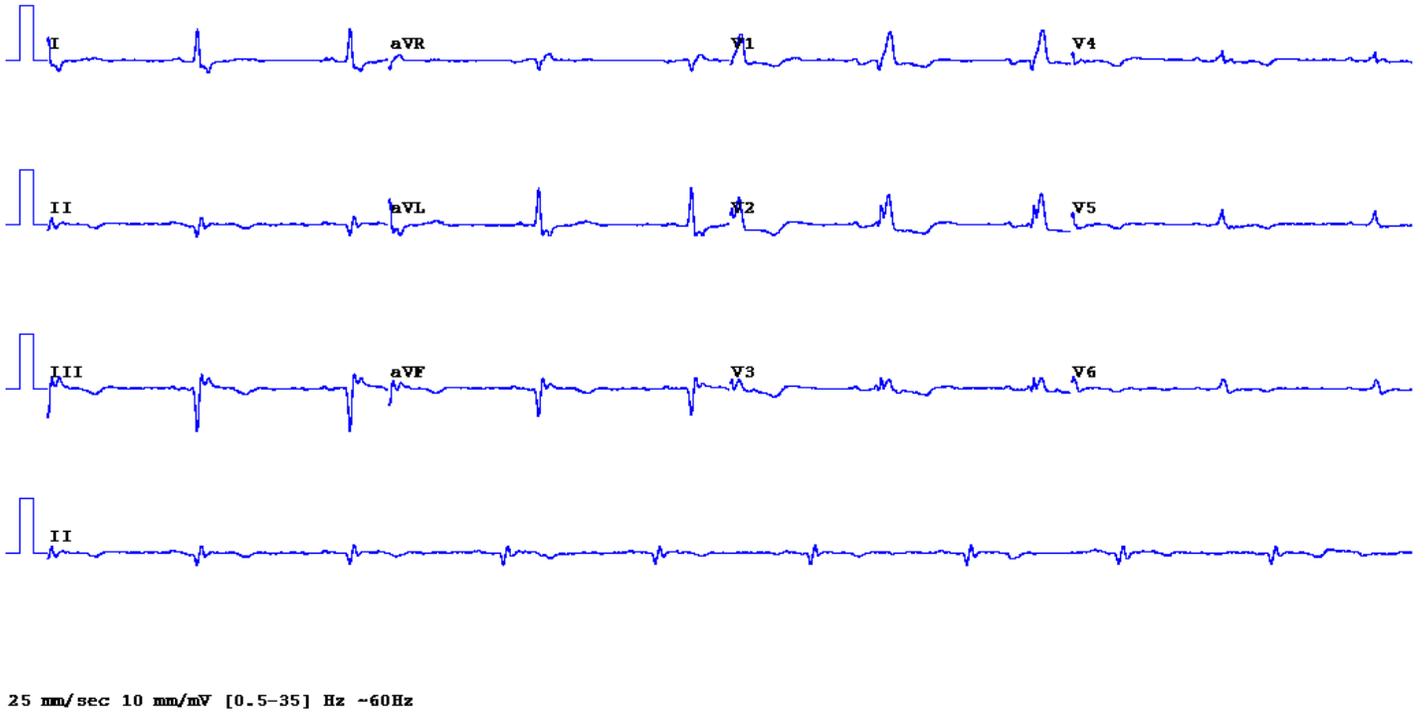


eFigure 3: Patient 02 old inferior MI STD ECG



eFigure 4: Patient 02 old inferior MI; the NEP ECG is similar.



eFigure 5: Patient 03 old inferior MI STD ECG

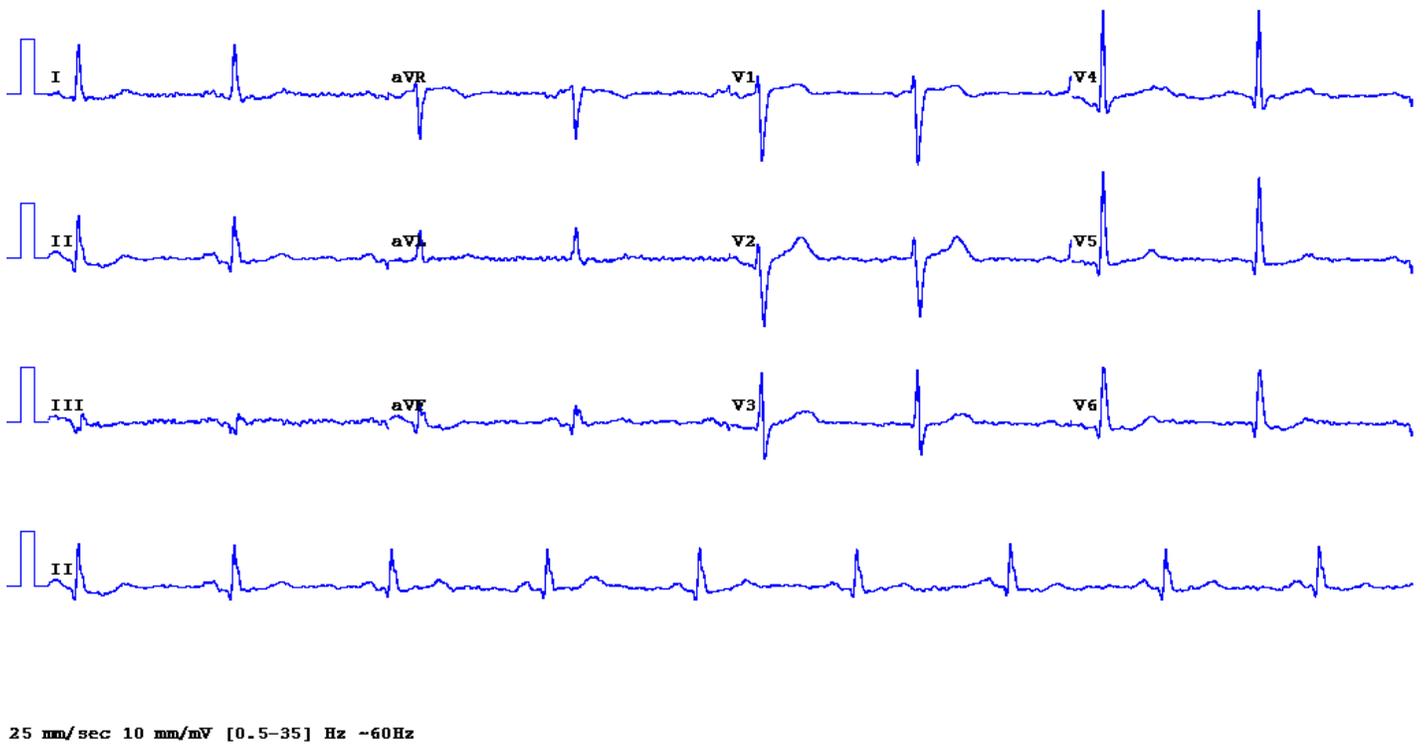

eFigure 6: Patient 03 old inferior MI; the NEP ECG is similar to the standard above.




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 7: Patient 04 old inferior and anterior MI STD ECG



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz


eFigure 8: Patient 04 old inferior MI NEP ECG is identical to standard ECG above.

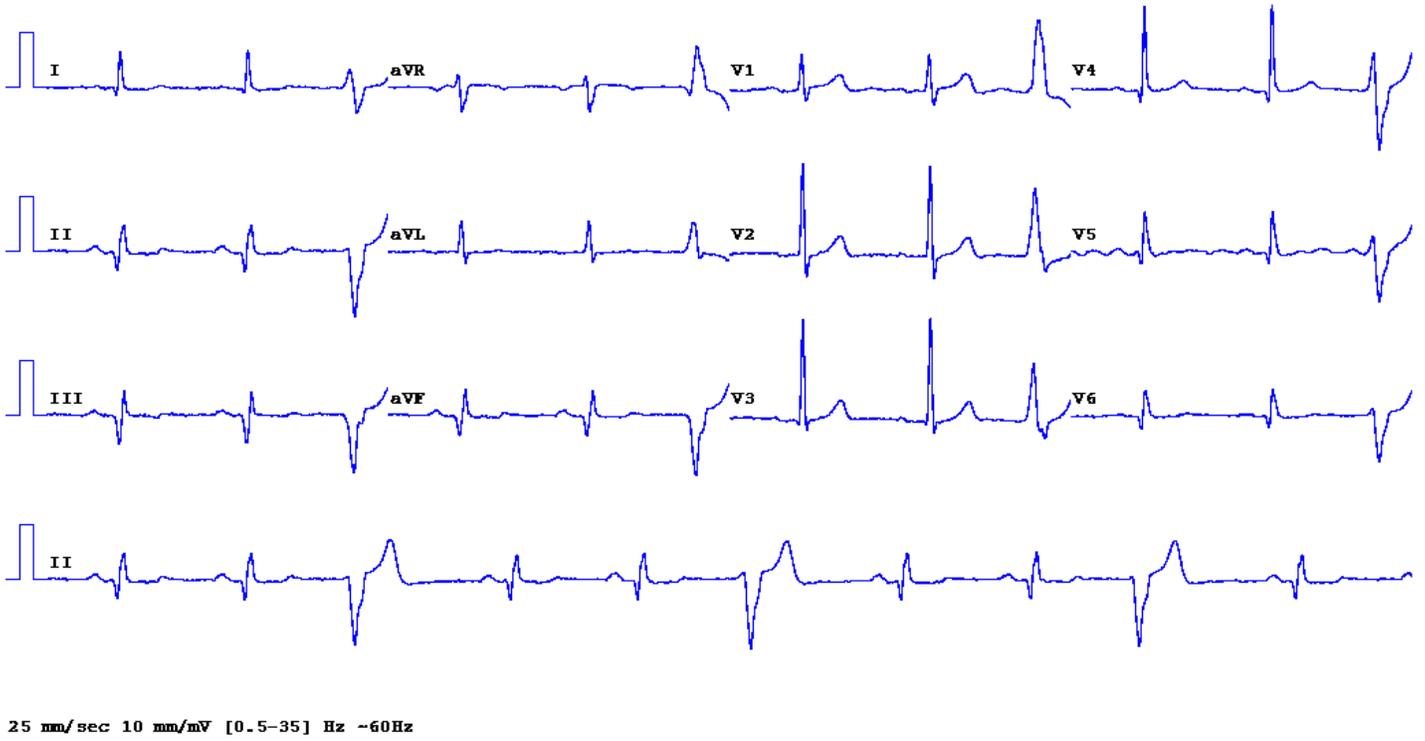

eFigure 9: Patient 05 old inferior MI STD ECG



eFigure 10: Patient 05 old inferior MI NEP ECG is similar.

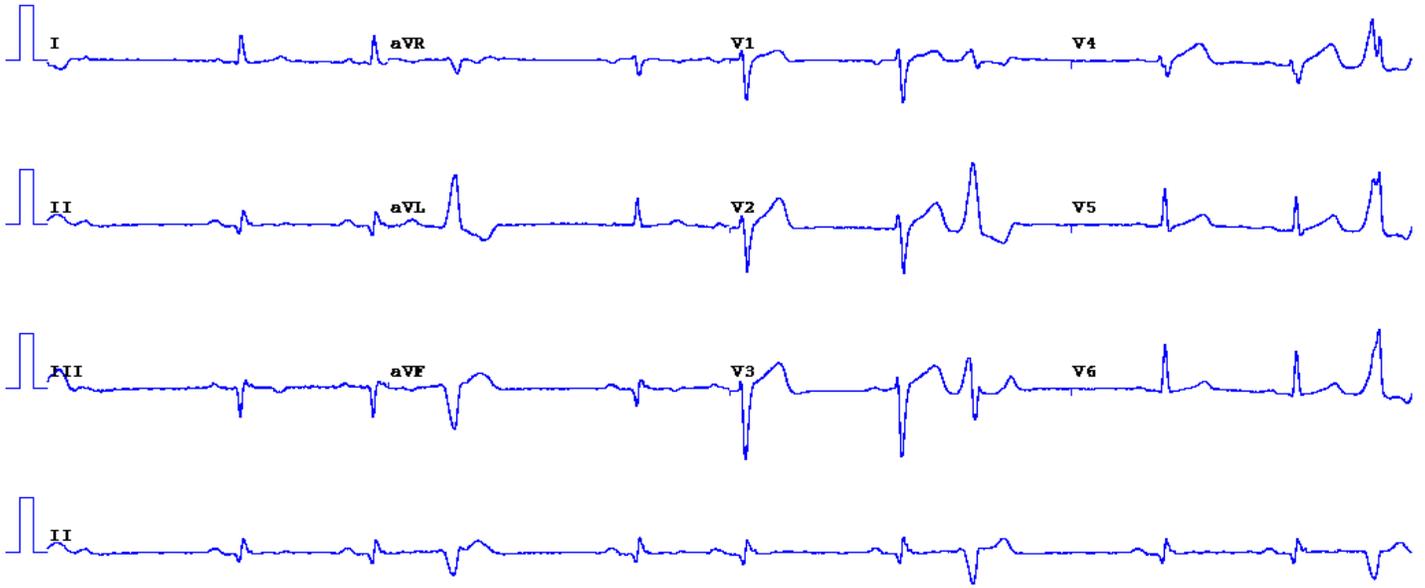

25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 11: Patient 06 old inferior MI, VPBs , STD ECG




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 12: Patient 06 old inferior MI; the NEP ECG is identical.




eFigure 13a: Patient 07 old inferior MI STD ECG.

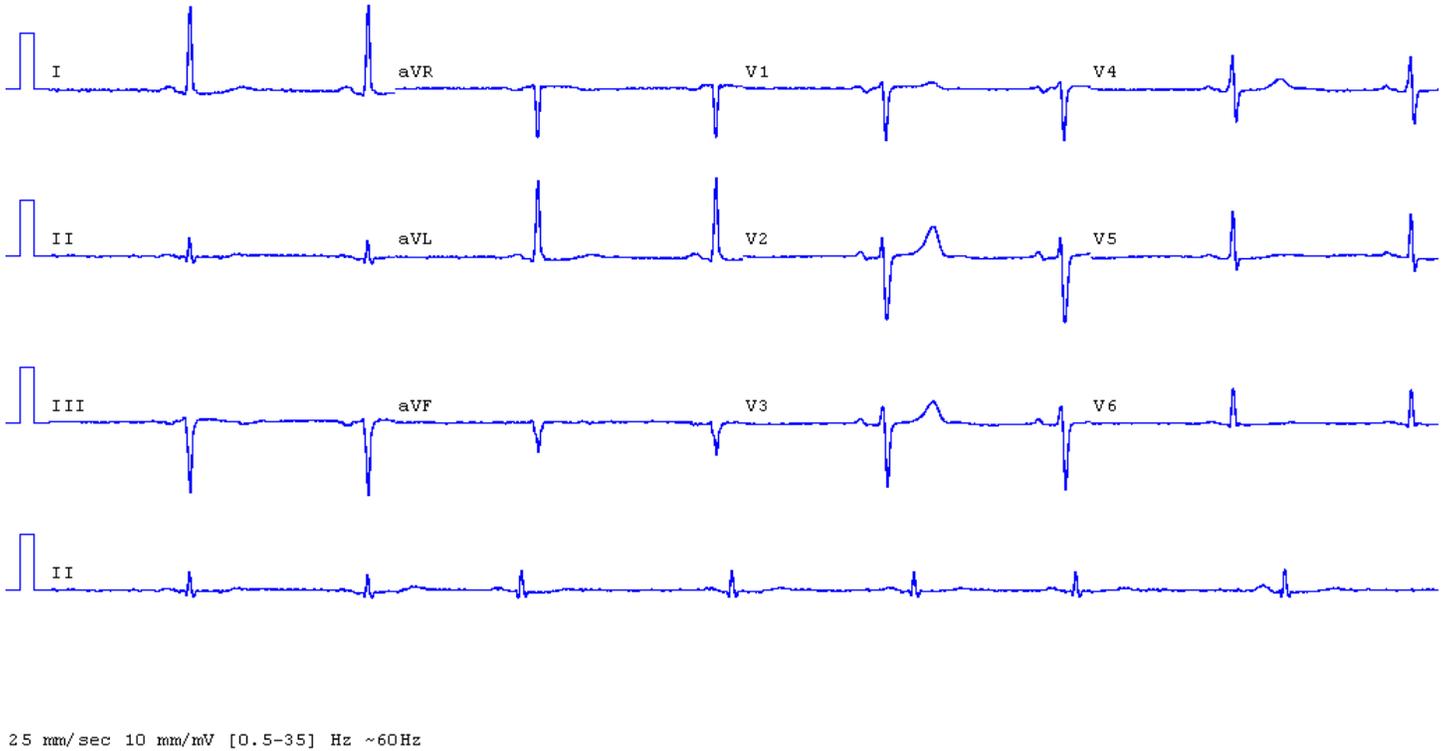

eFigure 13b: Patient 07 old inferior MI NEP ECG, similar.



eFigure 14: Patient 08 old inferior MI STD ECG

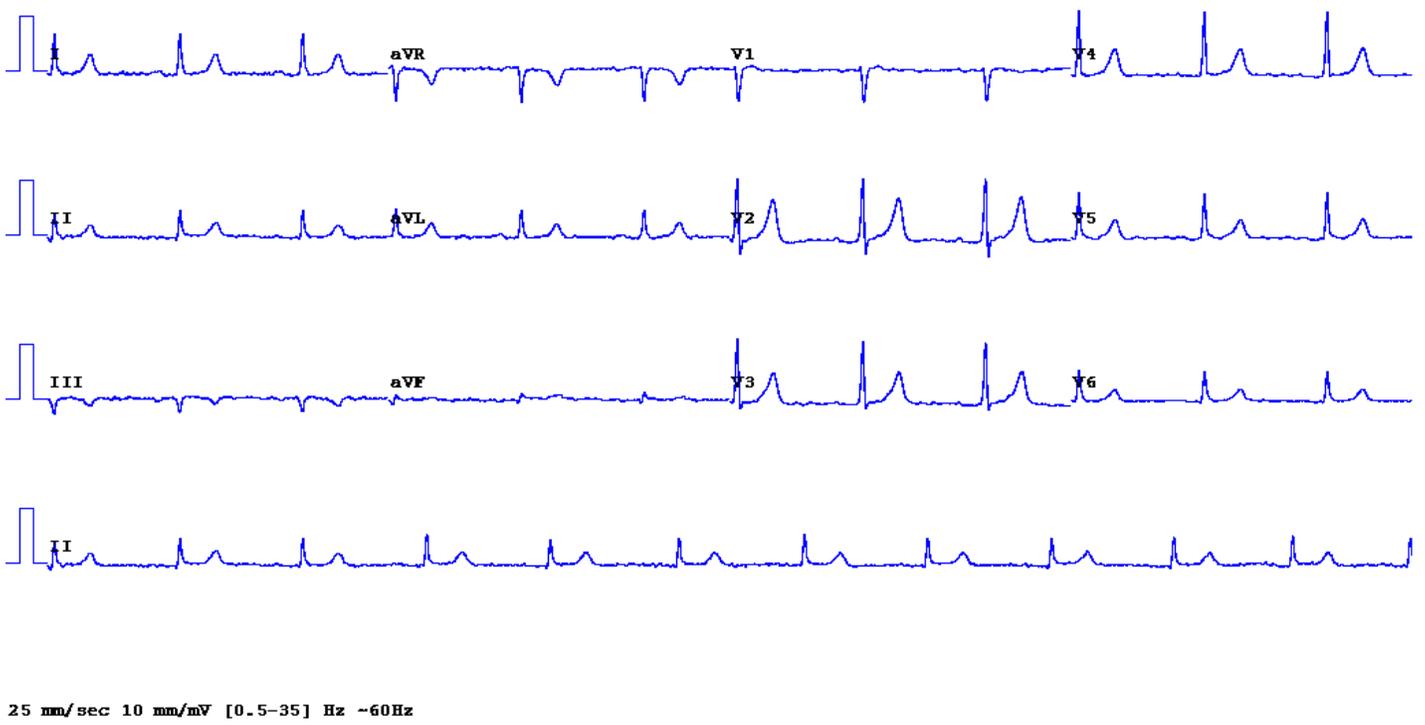

eFigure 15: Patient 08 old inferior MI NEP ECG is similar.




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 16: Patient 09 , inferior Q waves; old inferior MI STD ECG



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 17: Patient 09 old inferior MI , VPBs. NEP ECG is virtually identical.

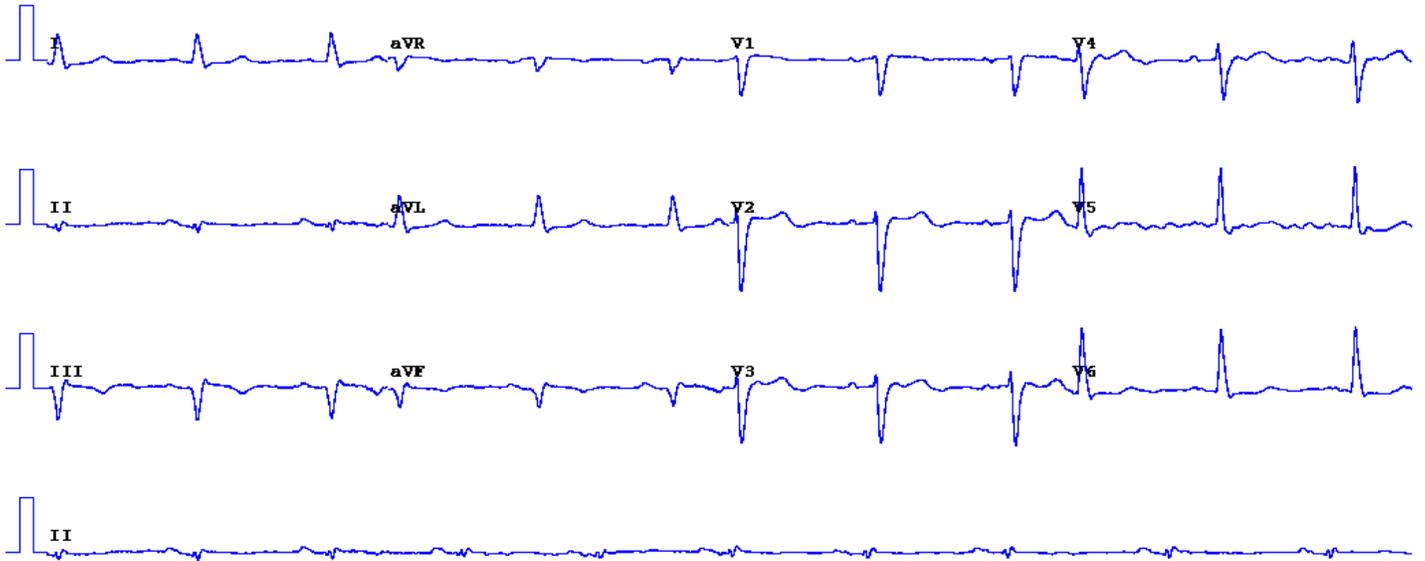

25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 18: Patient 10 old inferior MI STD ECG



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 19: Patient 10 old inferior MI NEP ECG is similar



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 20: Patient 11; deep inferior Q waves: old inferior MI STD ECG




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz


eFigure 21: Patient 11; similar Q waves and wave forms: old Inferior MI NEP ECG



eFigure 22: Patient 12 old inferior MI STD ECG, note artifacts, V4 – V 6. aVF.

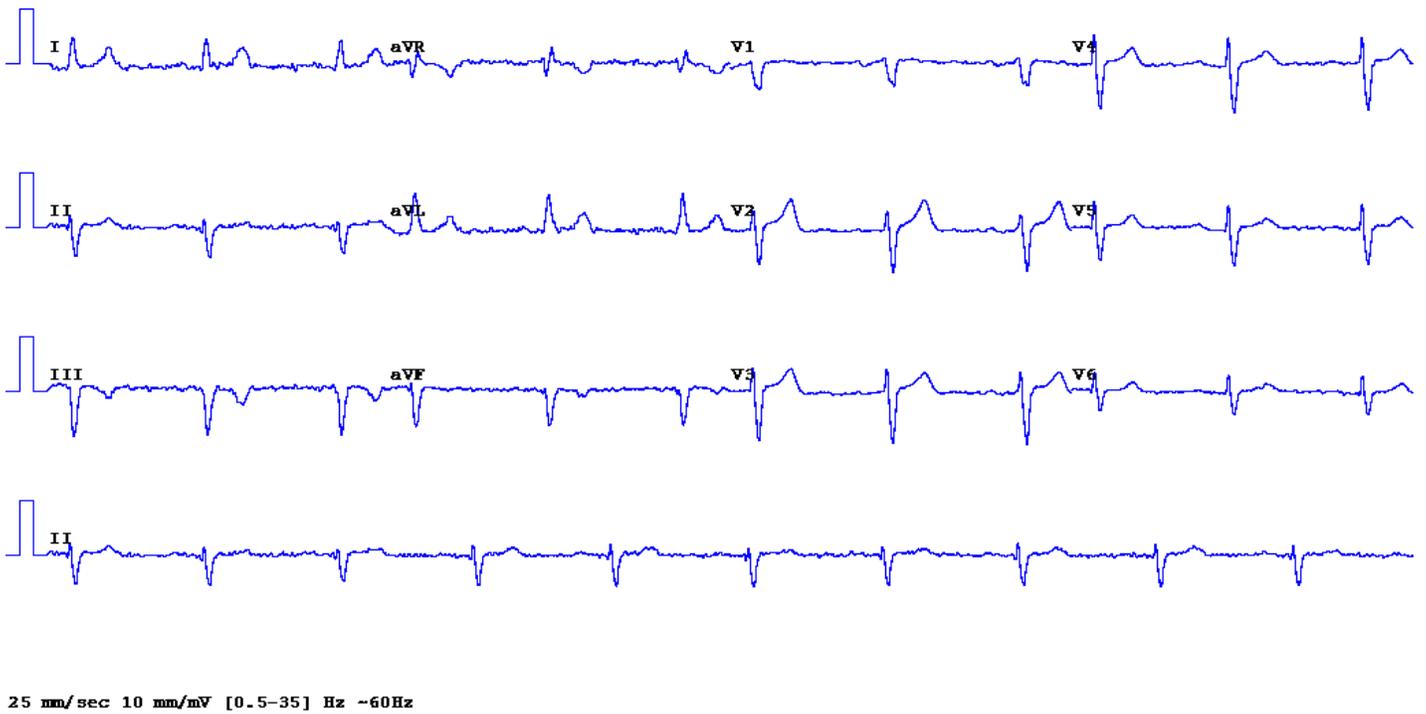

eFigure 23: Patient 12 old inferior MI NEP ECG shows clearing of artifacts. Better quality ECG.

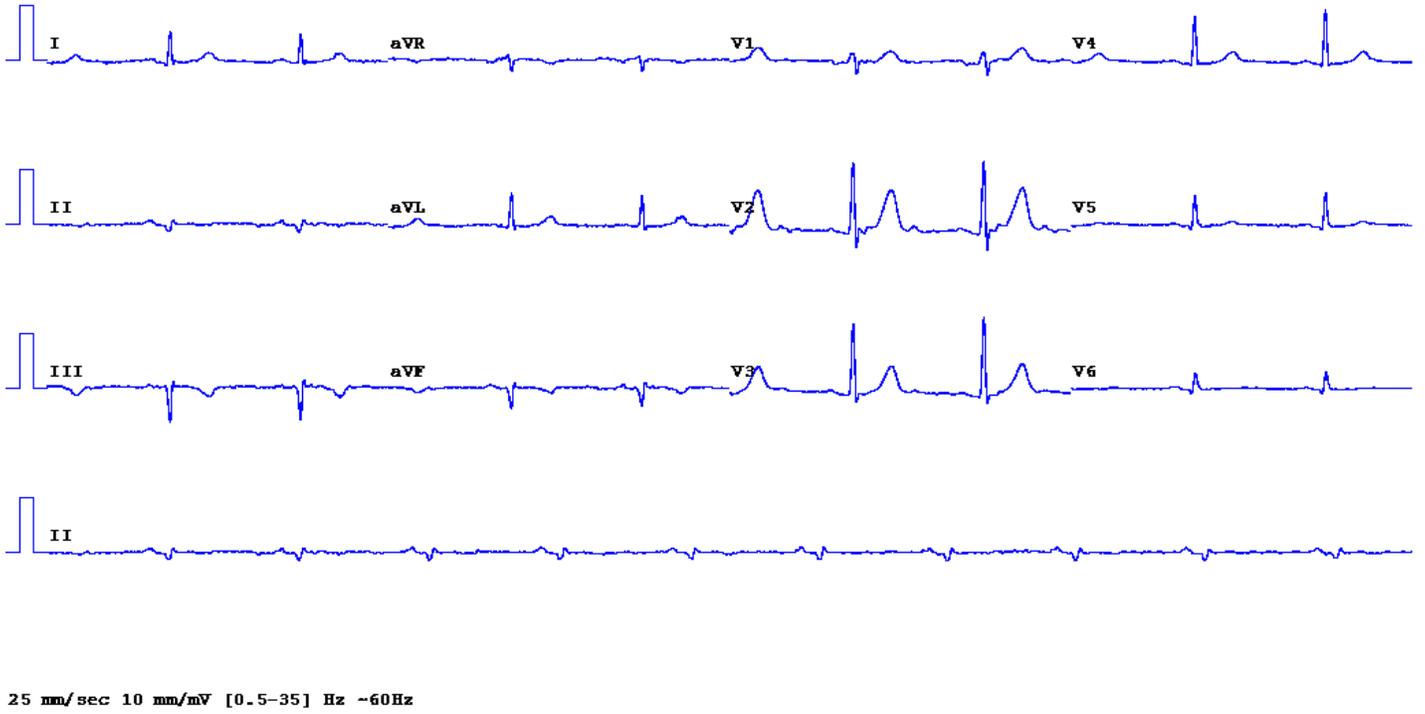

eFigure 24: Patient 13 old inferior MI STD ECG




eFigure 25: Patient 13 old inferior MI; the NEP ECG is virtually identical.




eFigure 26: Patient 14; deep wide Q waves,11,111, aVF: old inferior MI STD ECG




eFigure 27: Patient 14 old inferior MI;NEP ECG is not exempt from minor artifacts.

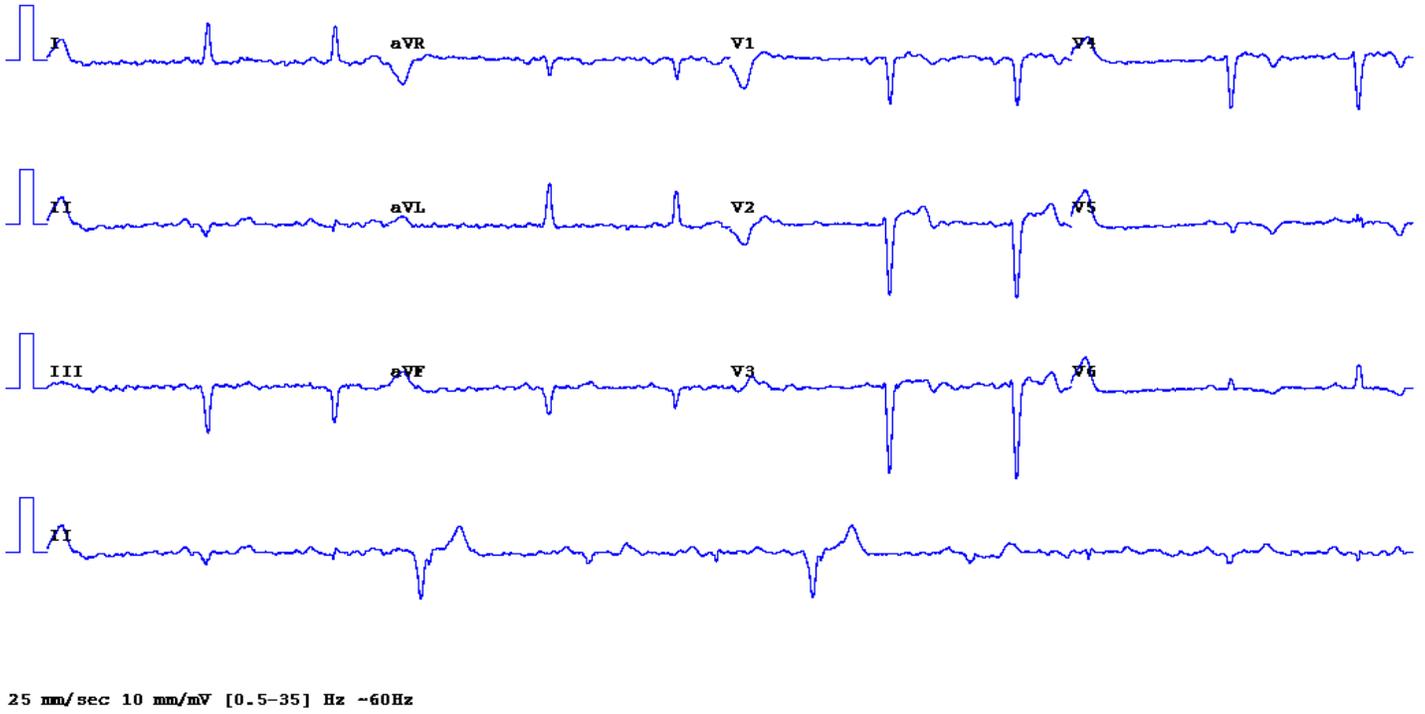


eFigure 28: Patient 15 old inferior MI STD ECG

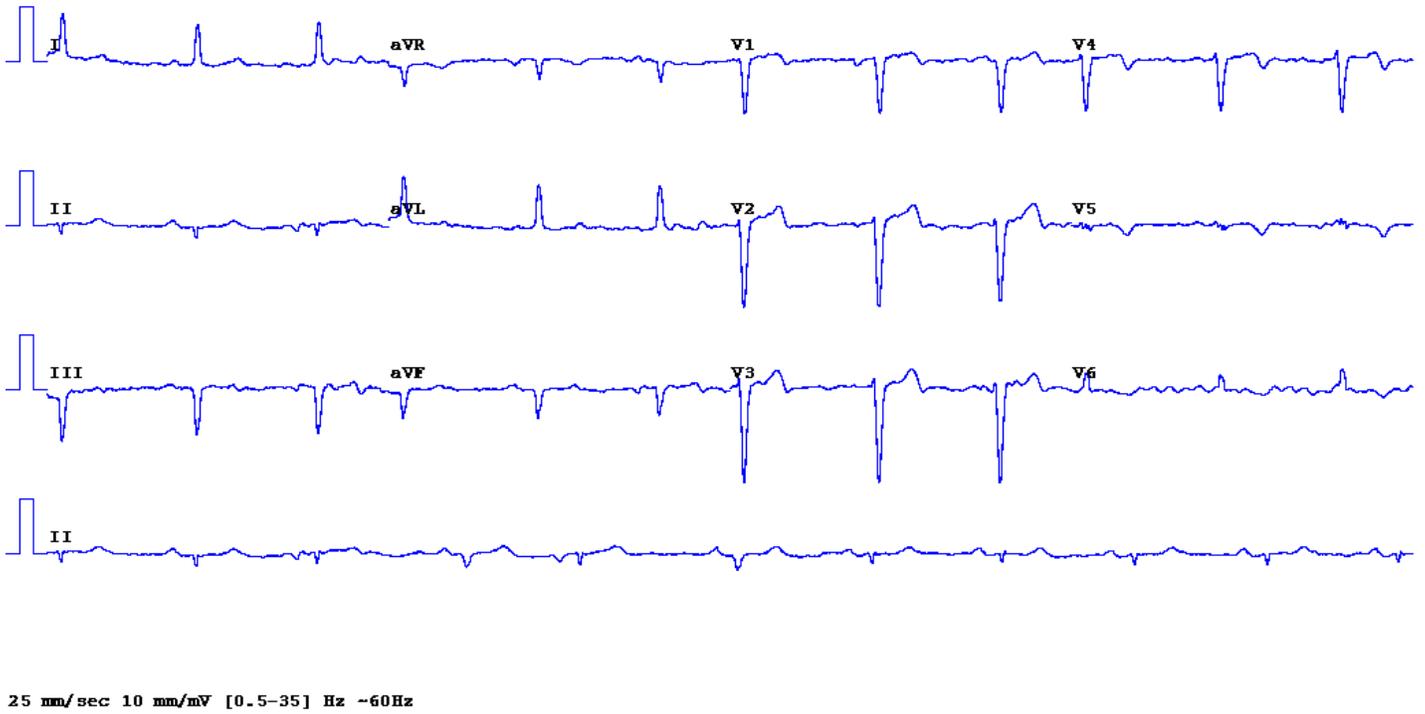


eFigure 29: Patient 15 old inferior MI; NEP ECG

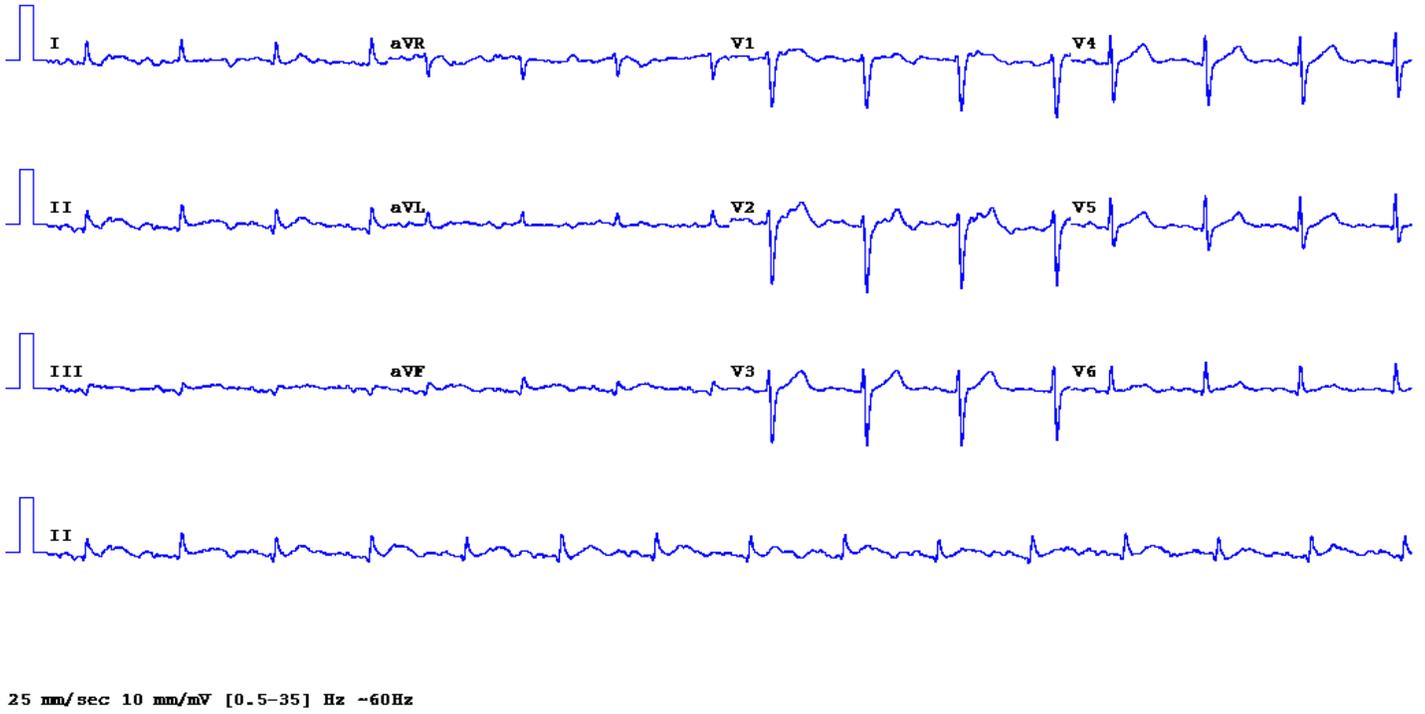


eFigure 30: Patient 16 old inferior MI STD ECG

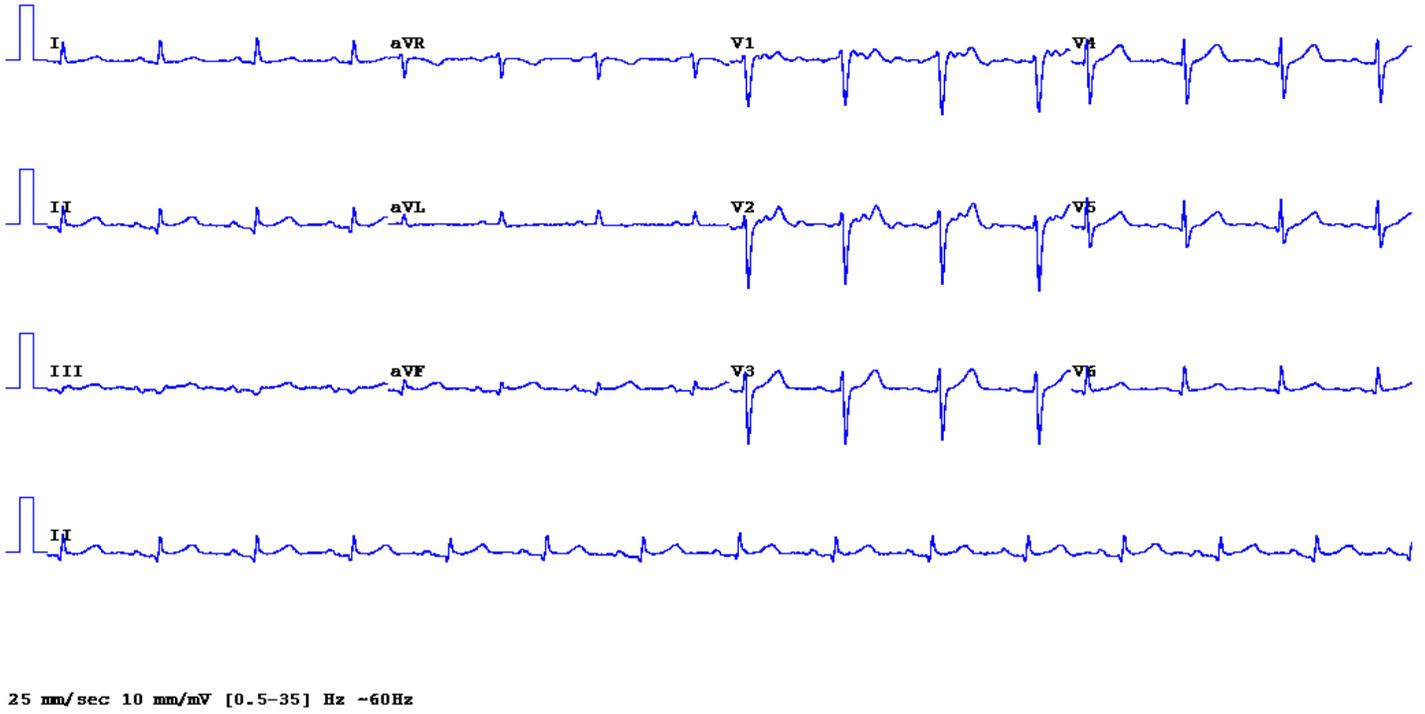

eFigure 31: Patient 16 old inferior MI; the NEP ECG is virtually identical.

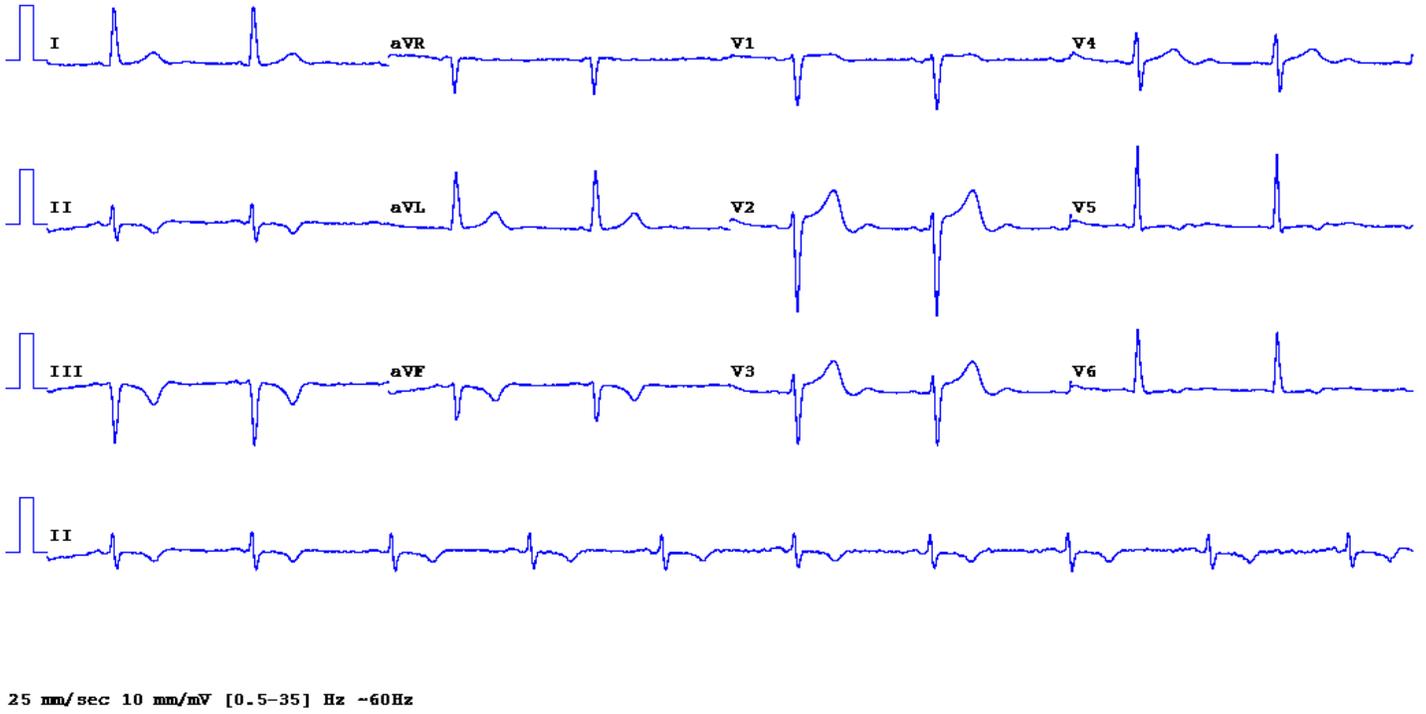

eFigure 32: Patient 17 old inferior MI NEP ECG

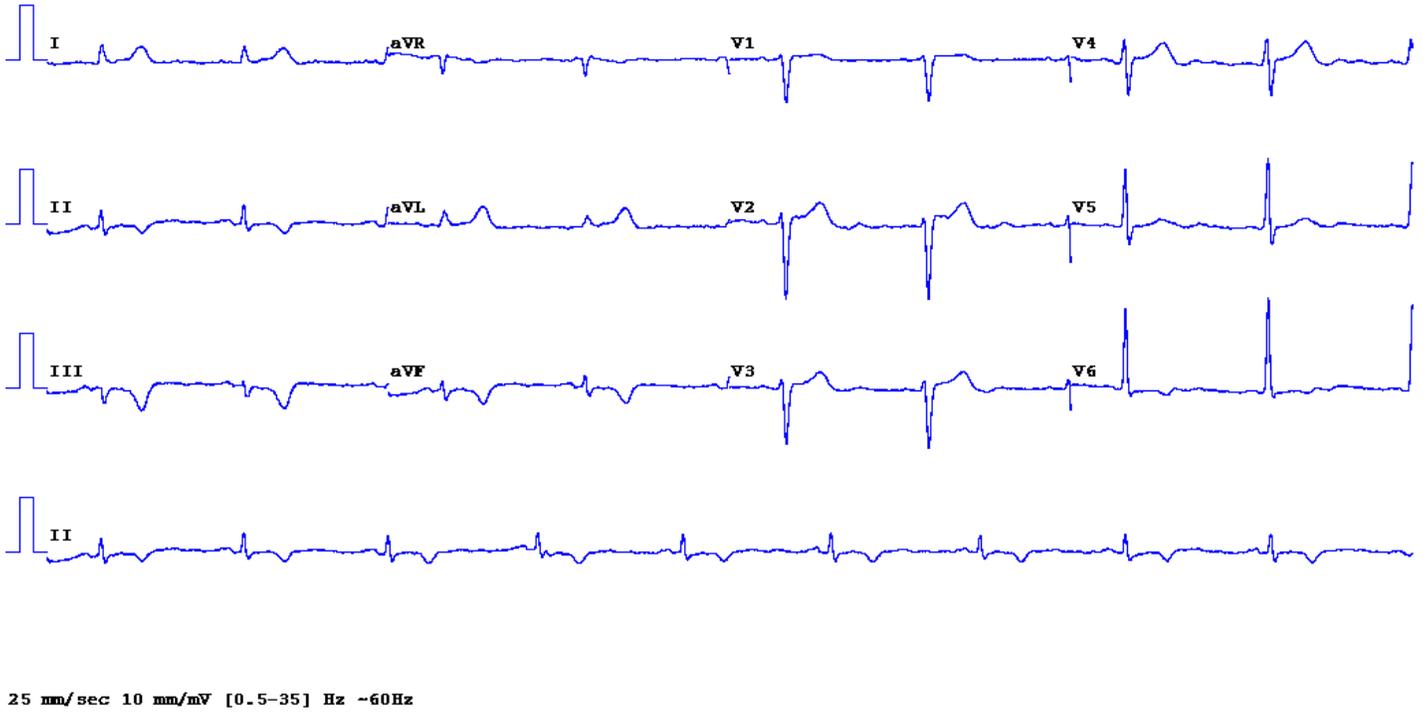


eFigure 33: Patient 17 old inferior MI; the STD ECG is similar.




eFigure 34: Patient 18 old inferior MI;STD ECG



eFigure 35: Patient 18 old inferior MI;NEP ECG

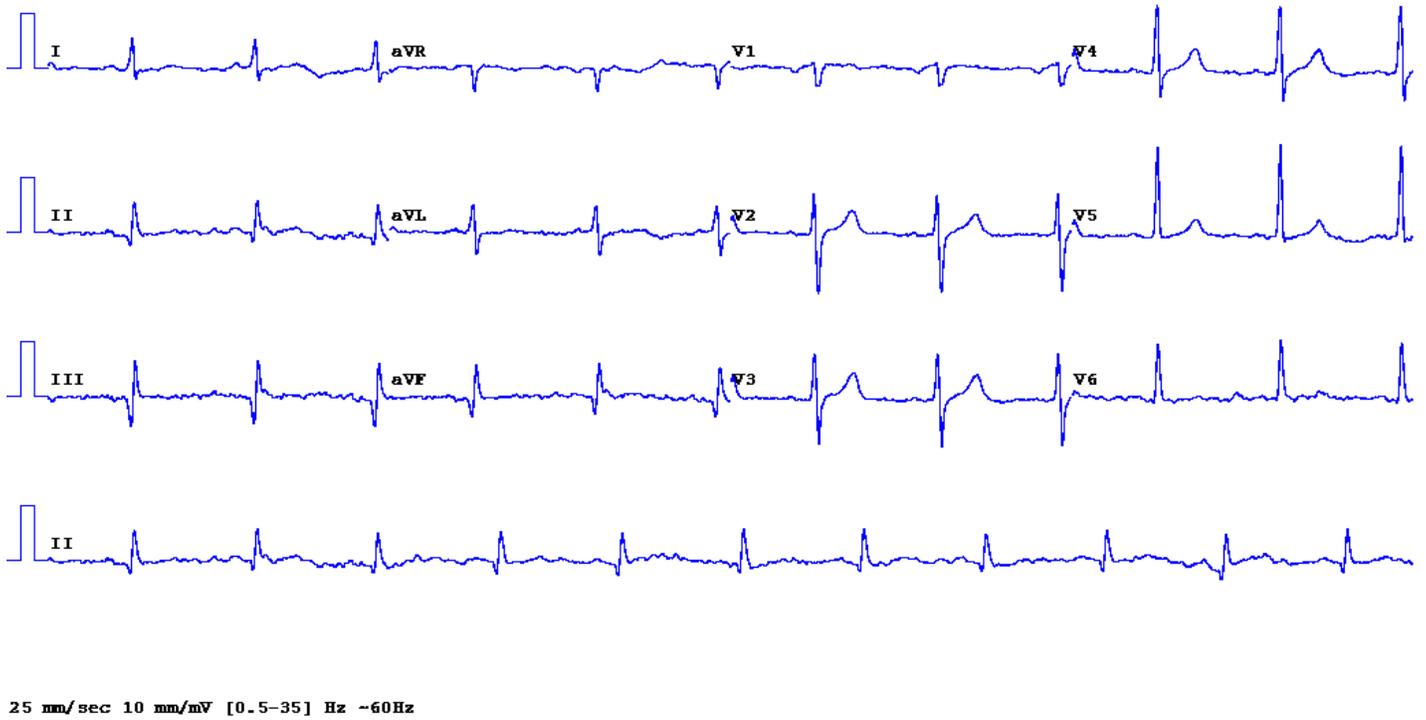


eFigure 36: Patient 19 old inferior MI:STD ECG

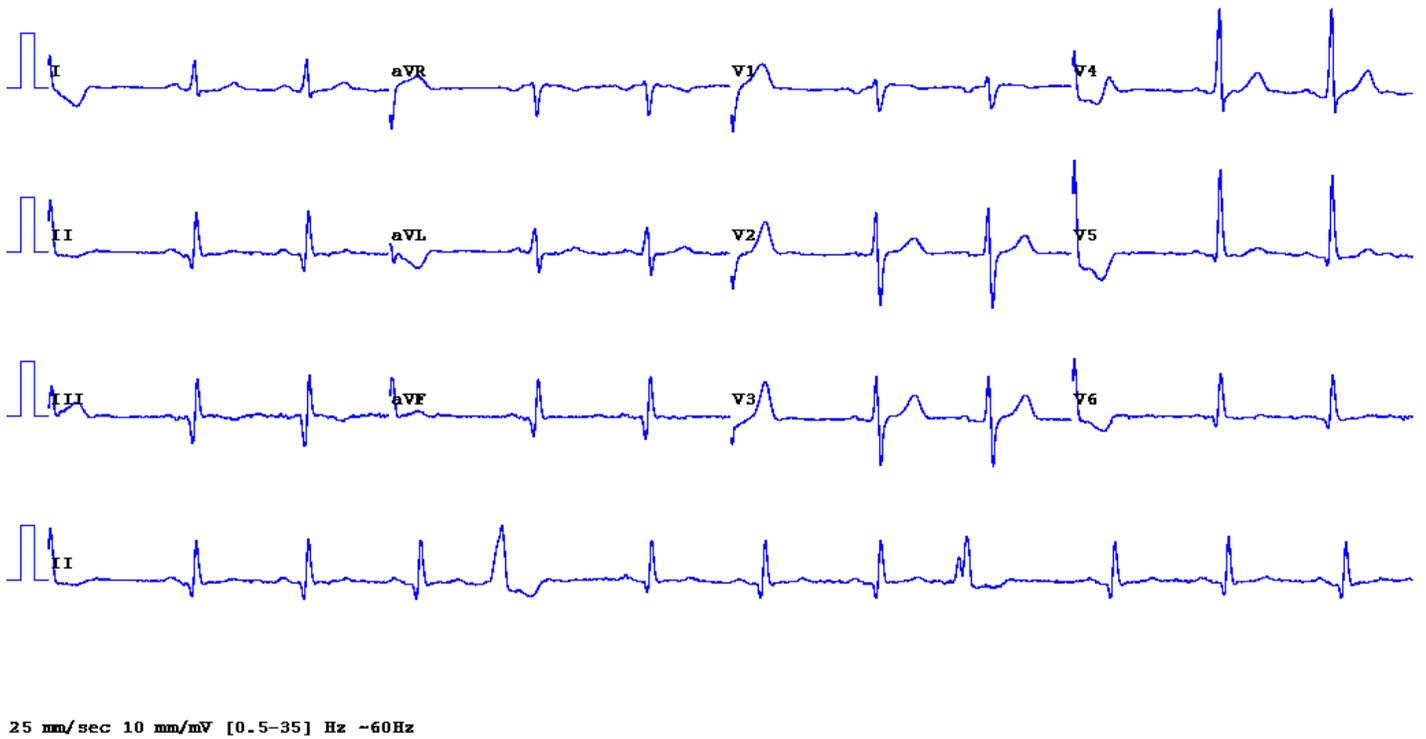

eFigure 37: Patient 19 old inferior MI: the NEP ECG is virtually identical.



eFigure 38: Patient 20 old inferior MI;STD ECG

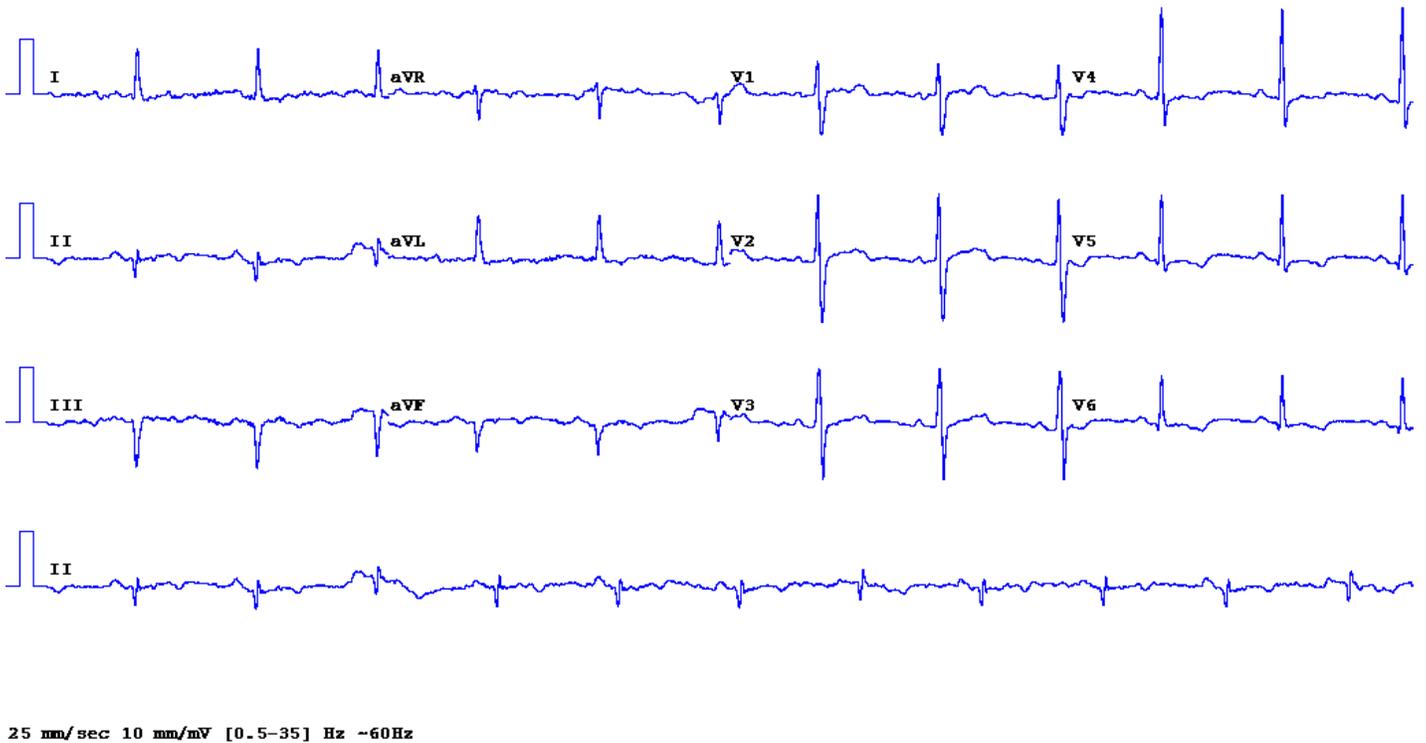

eFigure 39: Patient 20 old inferior MI: the NEP ECG is similar.



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 40: Patient 21 old inferior MI STD ECG




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 41: Patient 21 old inferior MI: the NEP ECG is similar.



eFigure 42: Patient 22 old inferior MI STD ECG

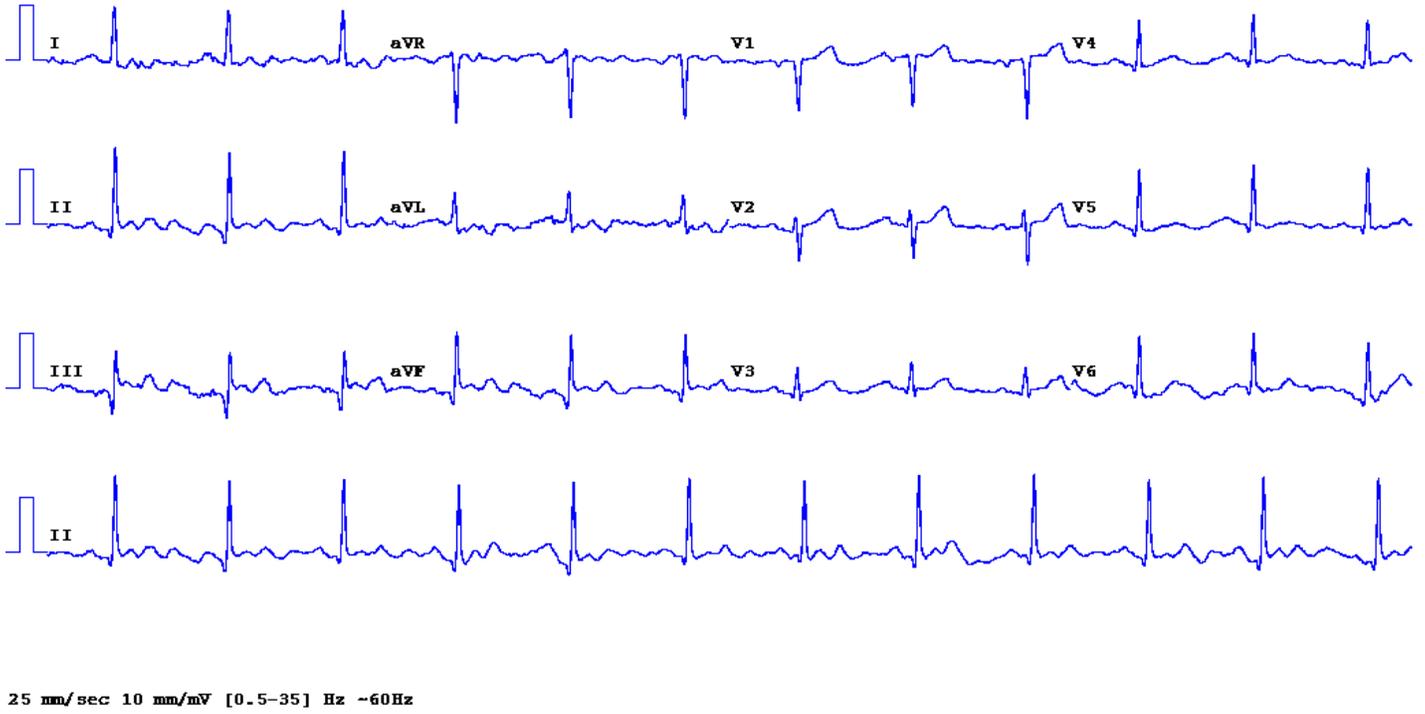

eFigure 43: Patient 22 old inferior MI the NEP ECG is similar.



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 44: Patient 23 old inferior MI STD ECG



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 45: Patient 23 old inferior MI; the NEP ECG is virtually identical.



eFigure 46: Patient 24 old inferior MI STD ECG



eFigure 47: Patient 24 old inferior MI; the NEP ECG is similar




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 48: Patient 25 inferior MI NEP ECG similar to STD




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz


eFigure 49: Patient 25 proven inferior MI STD ECG

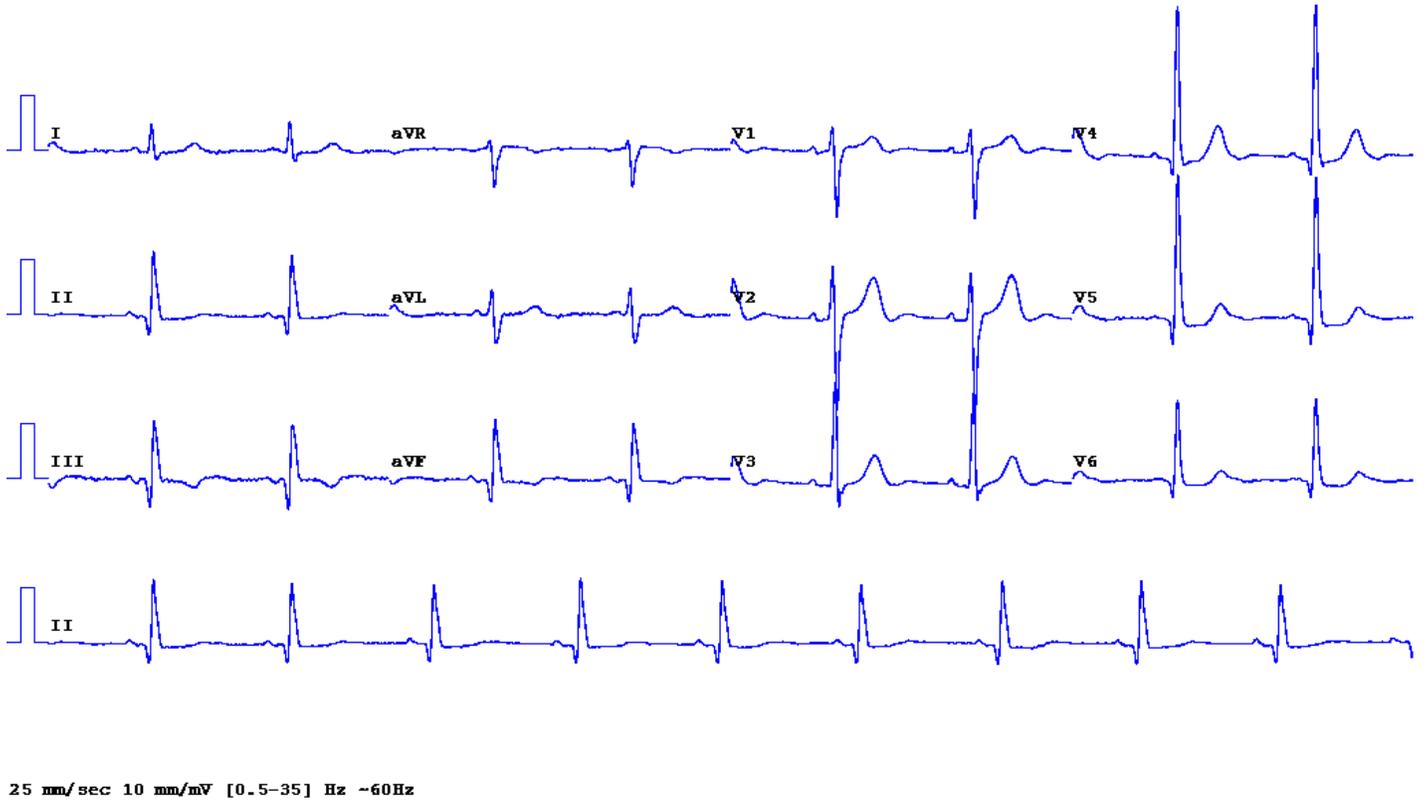

eFigure 50: Patient 26 inferior MI NEP ECG similar to STD

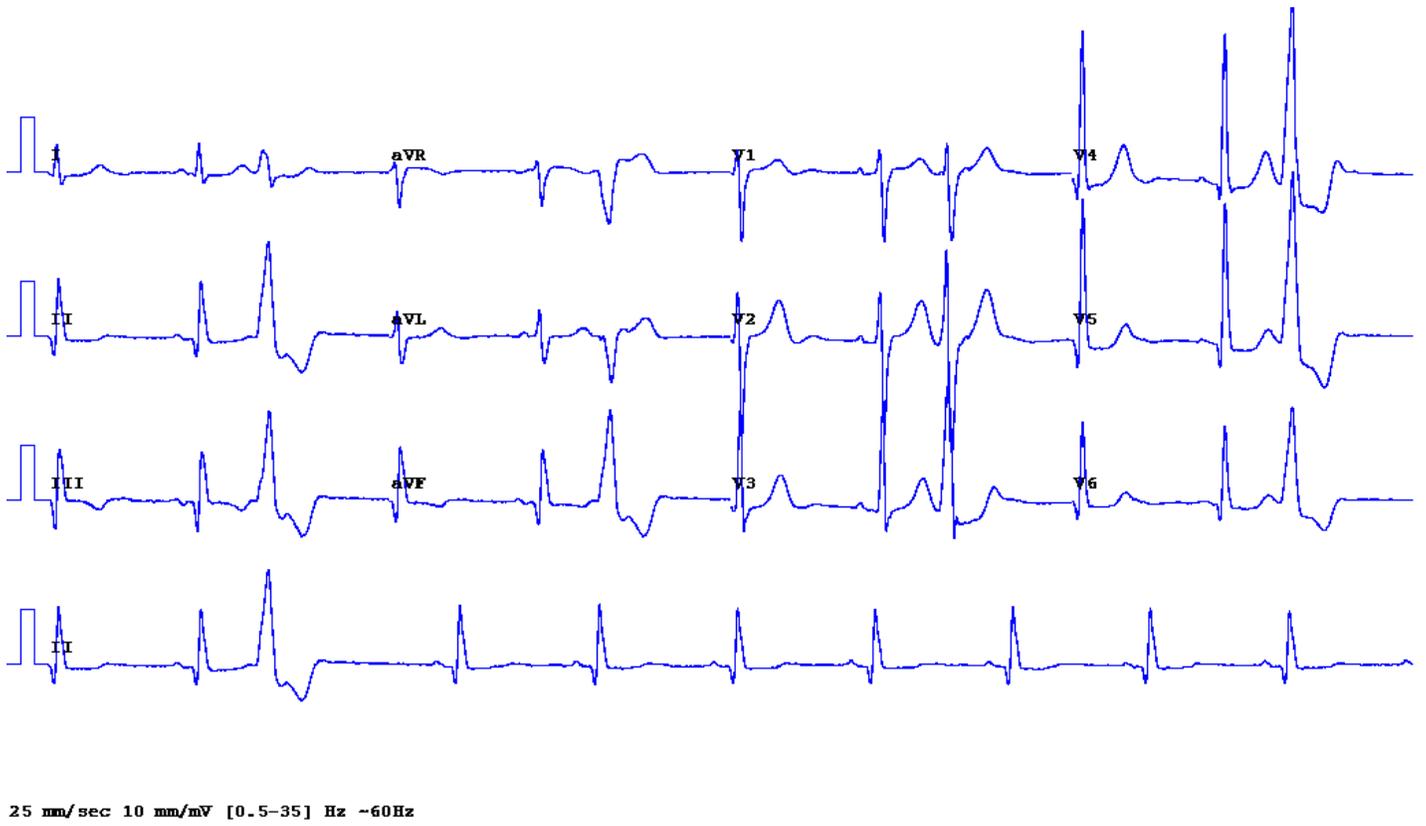


eFigure 51: Patient 26 proven inferior MI STD ECG

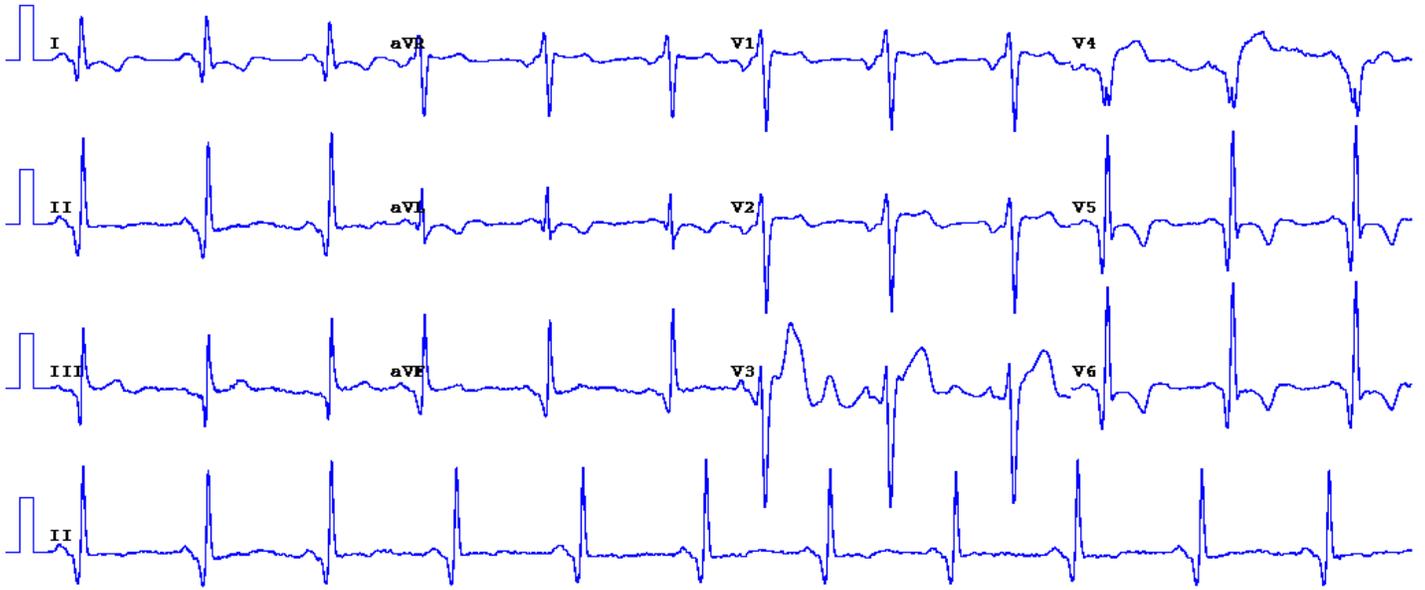

25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 52: Patient 27 old inferior MI STD ECG, note artifacts.

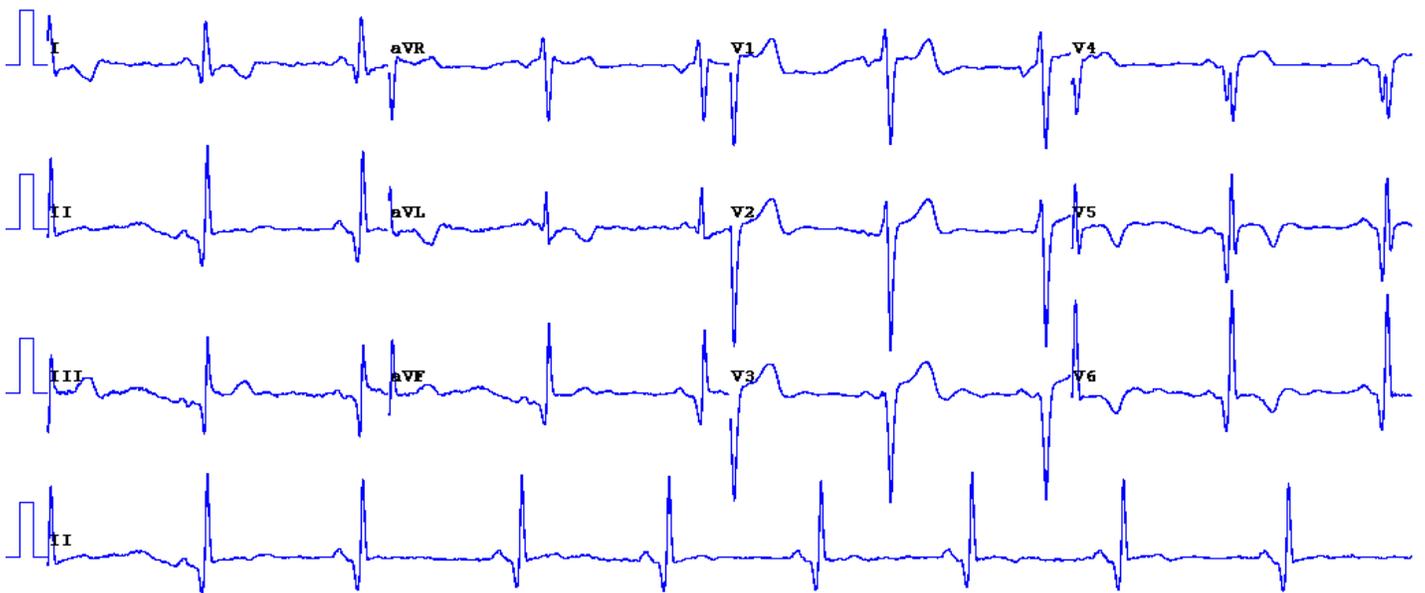

25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 53: Patient 27 old inferior MI NEP ECG , clearing of artifacts.



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 54: Patient 28 proven inferior MI STD ECG




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz


eFigure 55: Patient 28 inferior MI NEP ECG similar to STD




eFigure 56: Patient 29 old inferior MI STD ECG

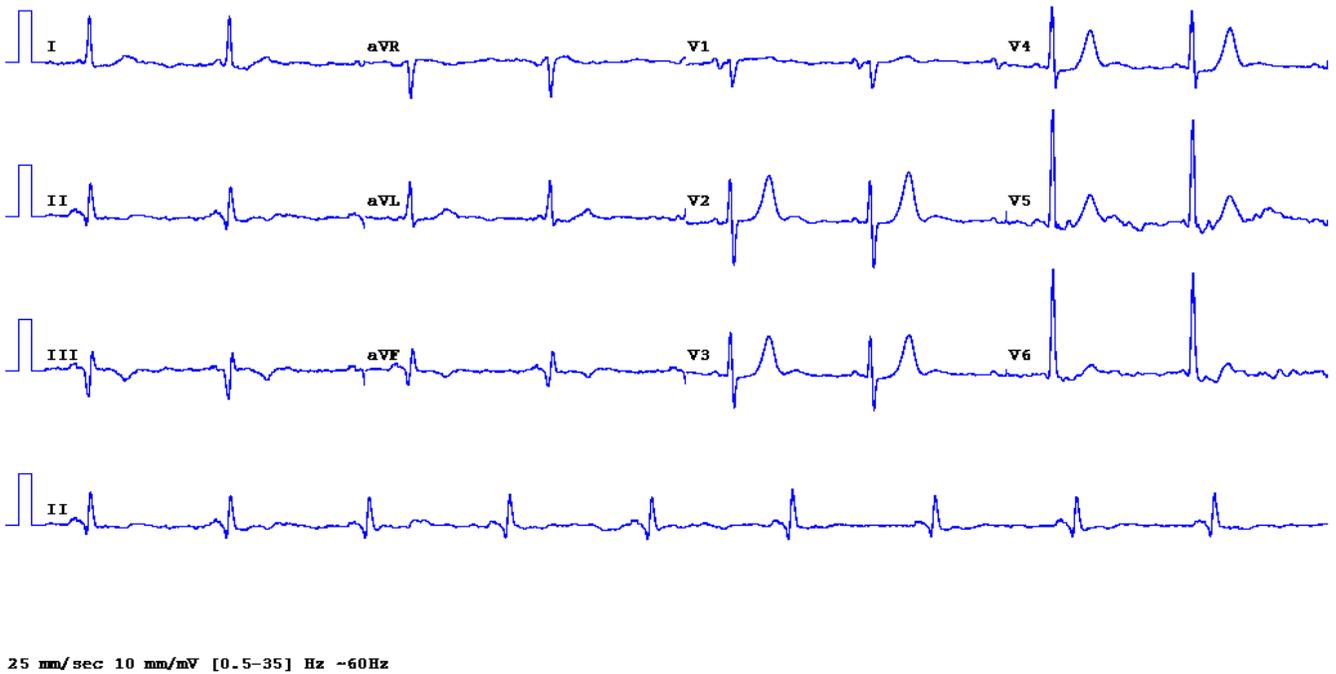

eFigure 57: Patient 29 old inferior MI; the NEP ECG is virtually identical

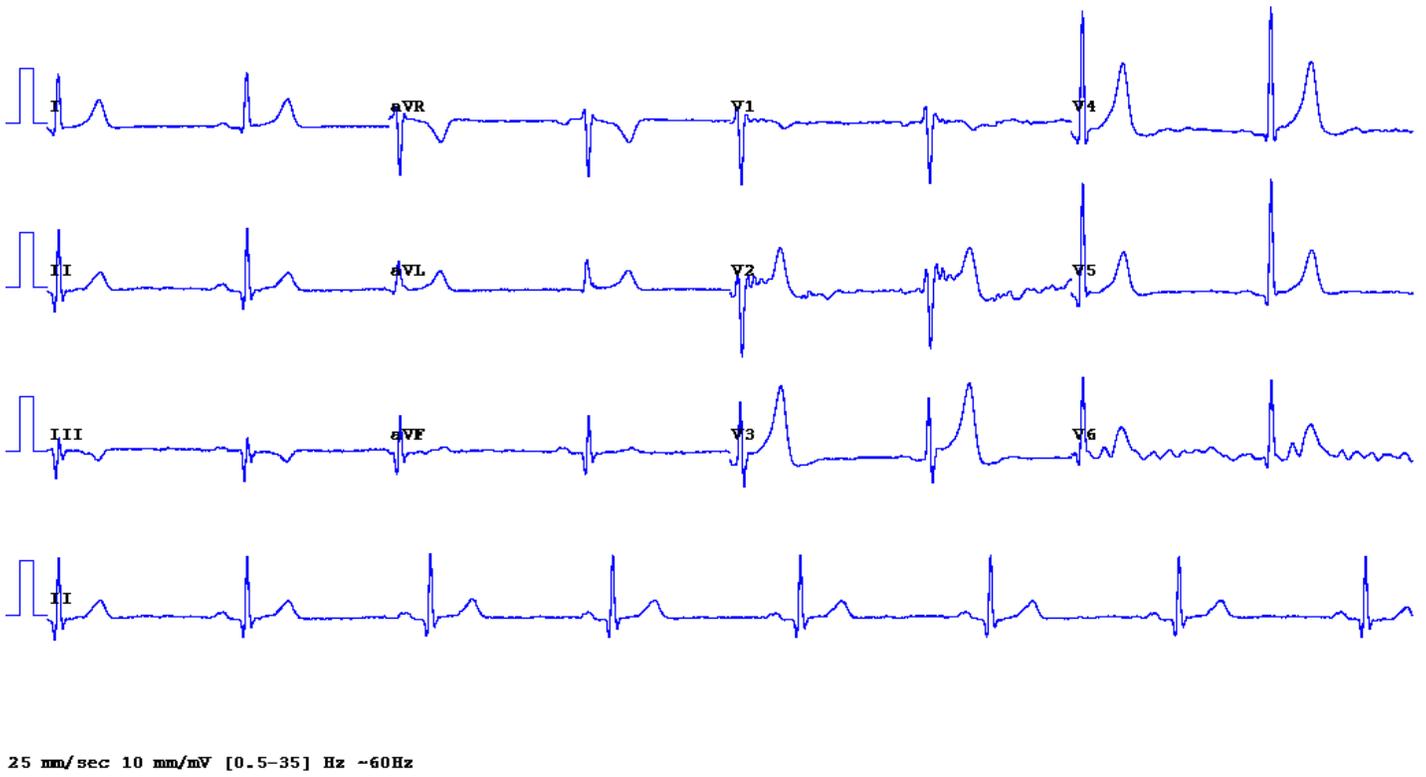

eFigure 58: Patient 30 old inferior MI STD ECG; artifacts

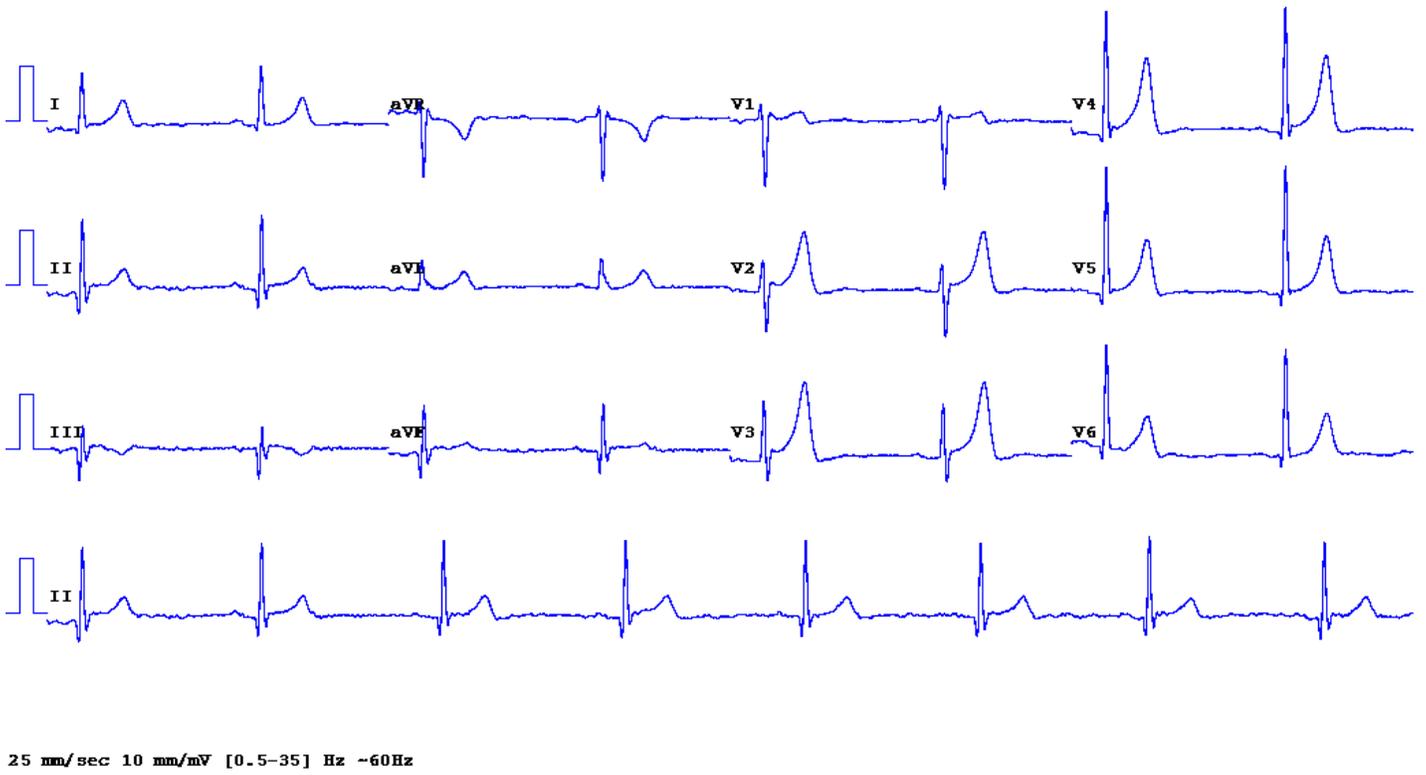


eFigure 59: Patient 30 old inferior MI NEP ECG artifacts V2 , V6 on standard ECG above cleared


## eAppendix 2




eFigure 1: Patient 01 low R wave amplitude, LEAD 1 AVL; STD ECG

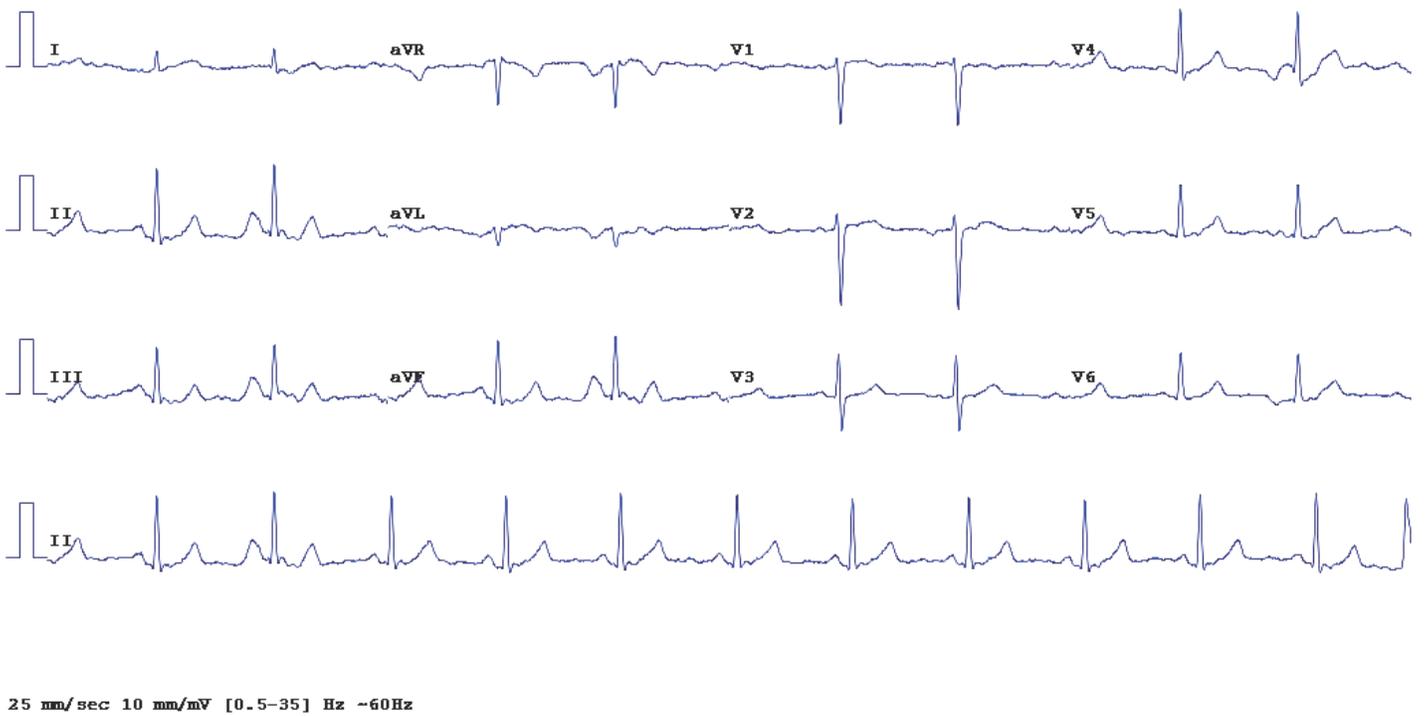


eFigure 2: Patient 01 LEAD 1 AVL NEP ECG R wave amplitude is similar to STD ECG

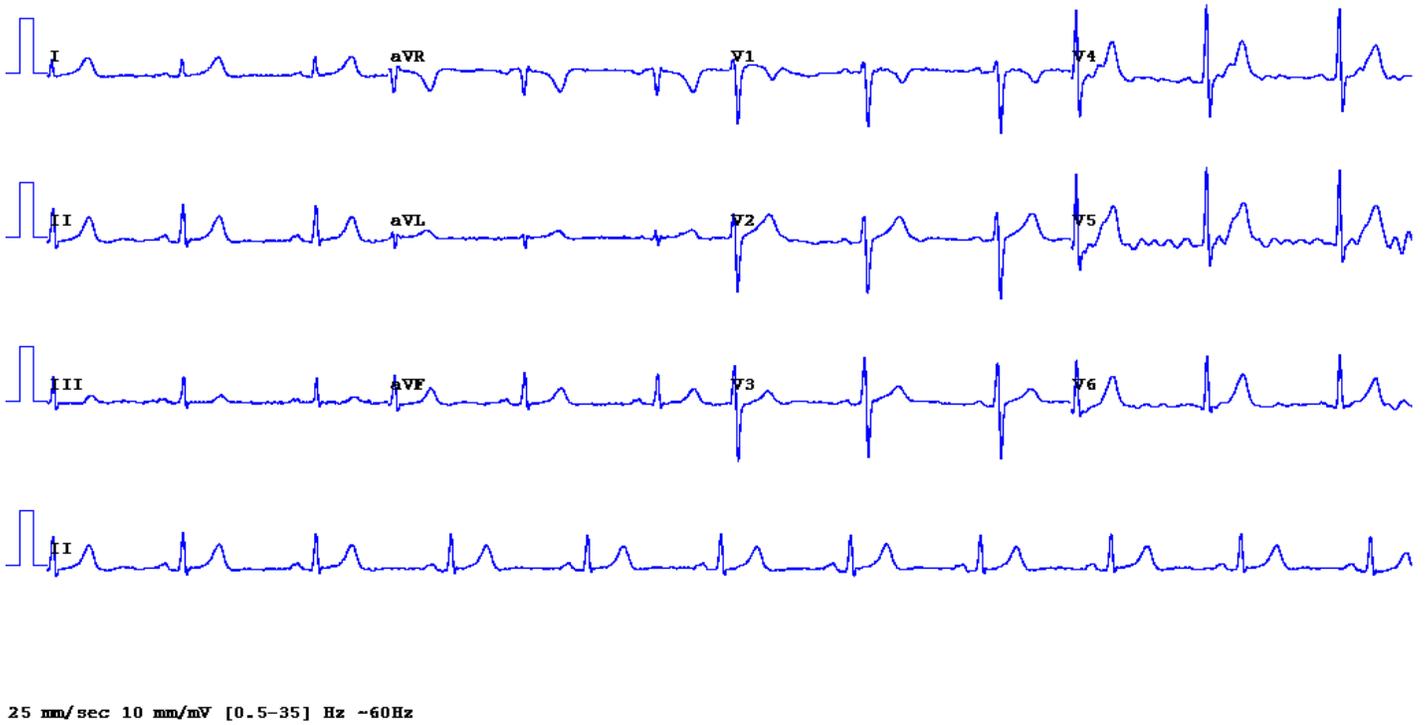

eFigure 3: Patient 02; Low R wave amplitude, lead 1 AVL; STD ECG.

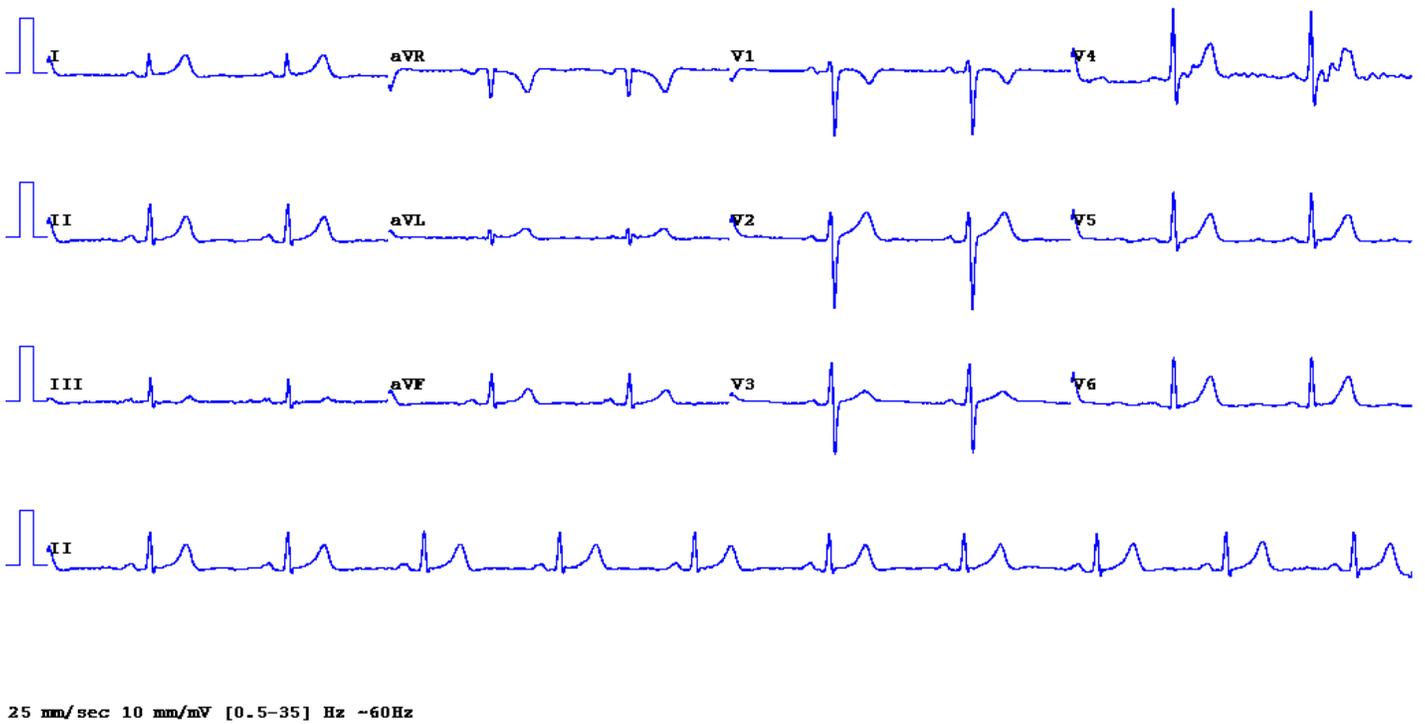

eFigure 4: Patient 02 LEAD 1 AVL NEP ECG, R wave amplitude in leads 1, aVL and other leads is similar to the STD ECG

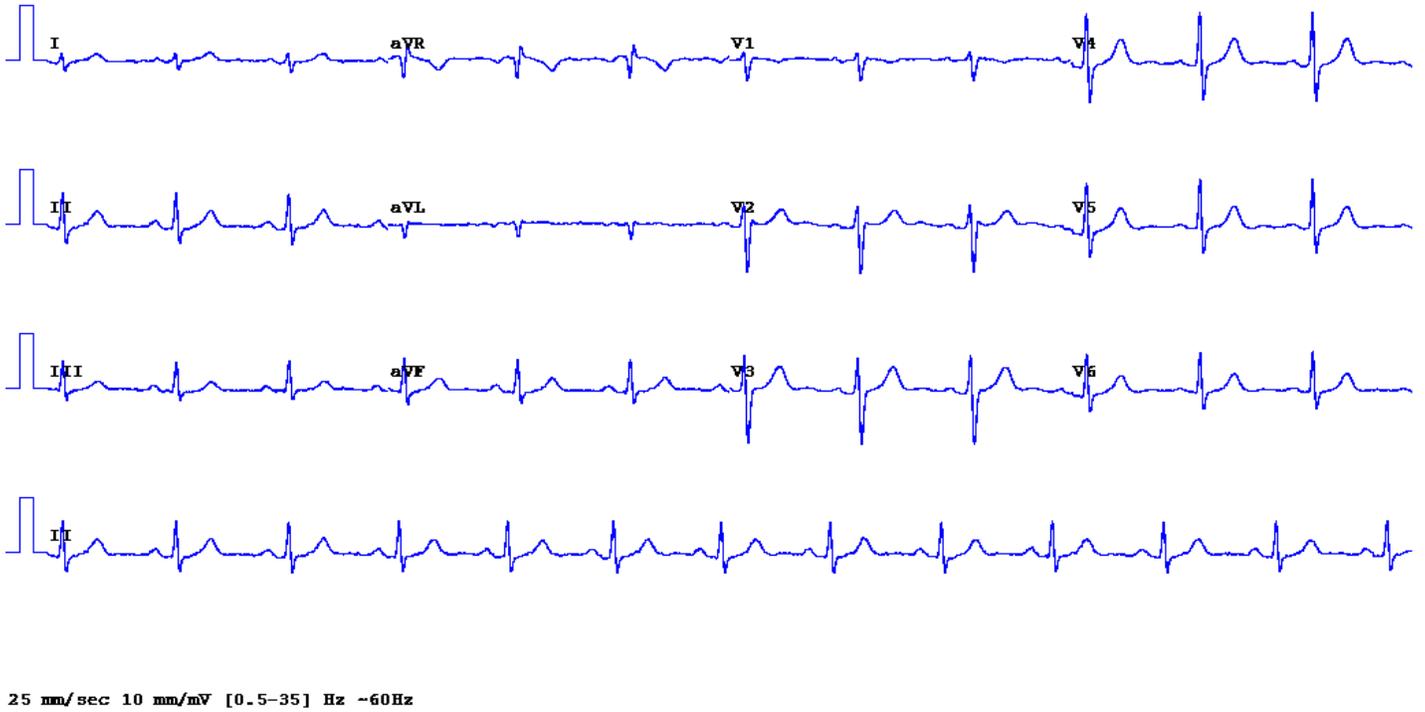

eFigure 5: Patient 03 STD ECG: very low R wave amplitude leads 1 avl



eFigure 6: Patient 03 Lead 1 AVL NEP ECG no change in very low R Wave amplitude compared with STD; no conversion to QS pattern.




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 7: Patient 04 LEAD 1 AVL STD ECG



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 8: Patient 04 LEAD 1 AVL NEP ECG similar

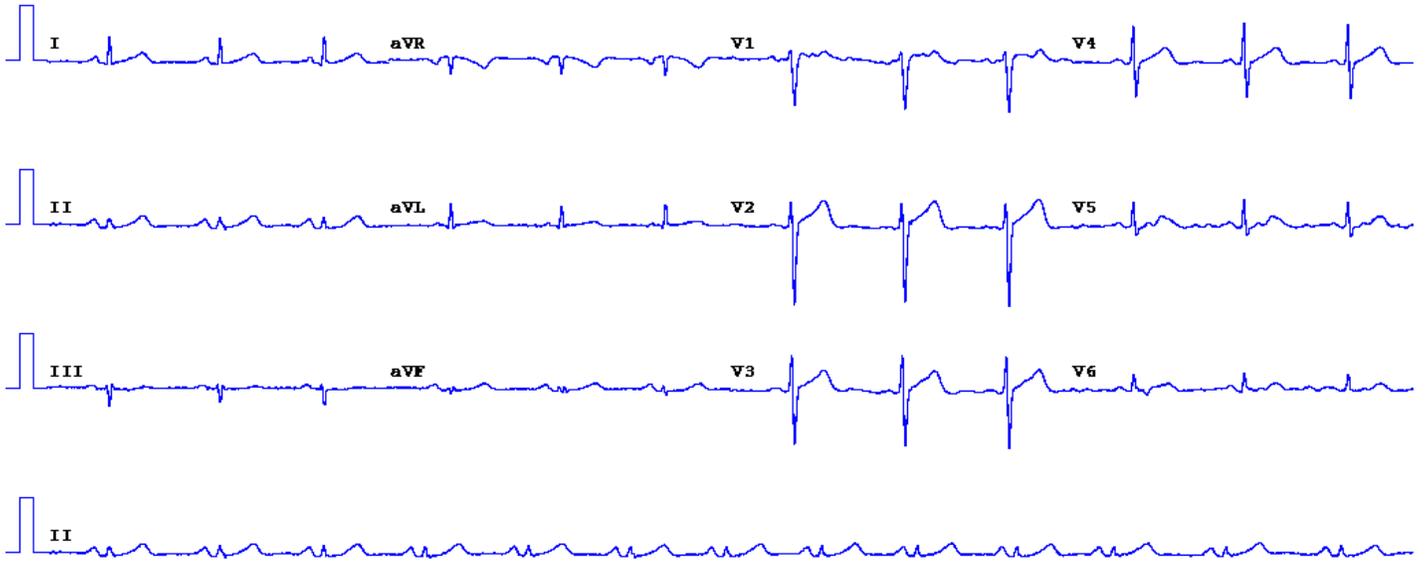

25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 9: Patient 05 LEAD 1 AVL STD ECG



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 10: Patient 05 LEAD 1 AVL NEP ECG R wave amplitude in all leads are similar to STD



eFigure 11: Patient 06 LEAD 1 AVL STD ECG; very low amplitude R waves 1 and aVL




eFigure 12: Patient 06 LEAD 1 AVL NEP ECG R wave amplitude in all 12 leads similar to STD



eFigure 13: Patient 07 LEAD 1 AVL NEP ECG all R waves similar to STD ECG




eFigure 14: Patient 07 LEAD 1 AVL STD ECG

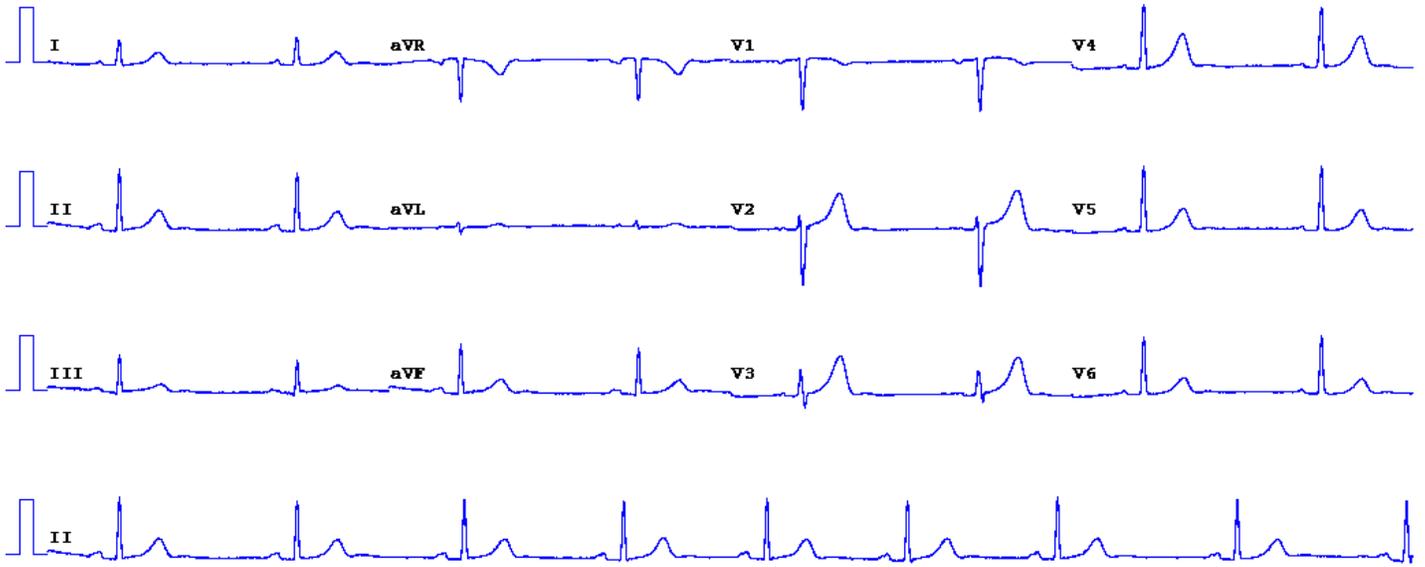


25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 15: Patient 08 LEAD 1 AVL STD ECG

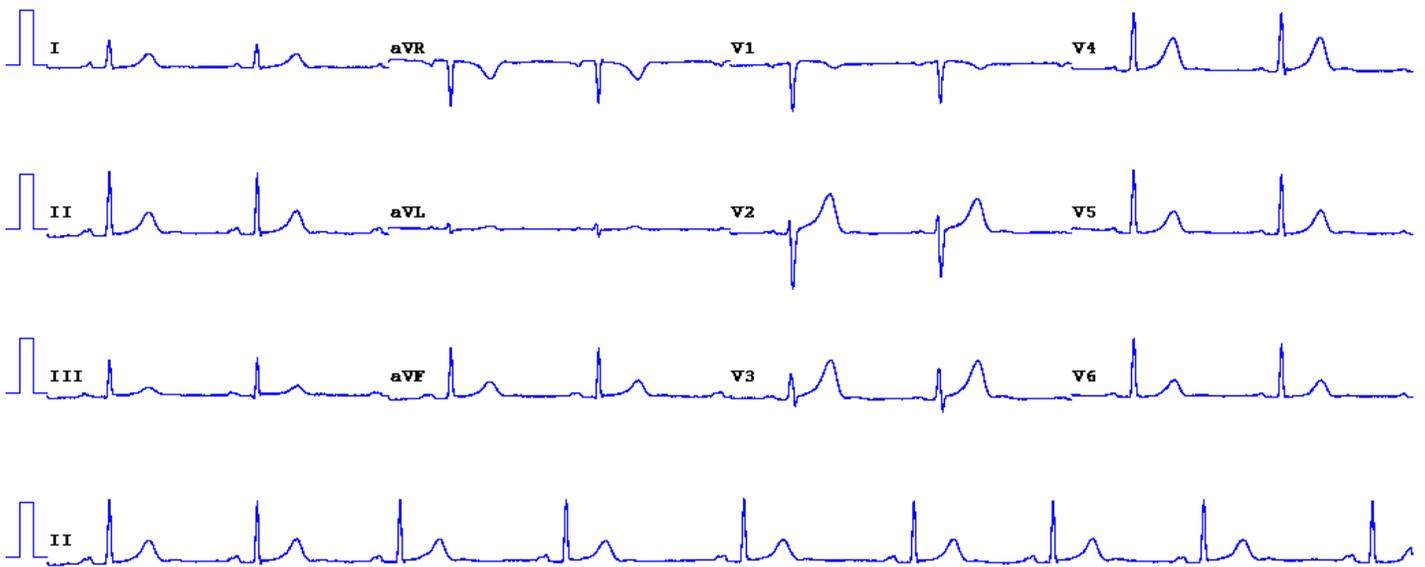

25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 16: Patient 08 LEAD 1 AVL NEP ECG , R waves similar to STD

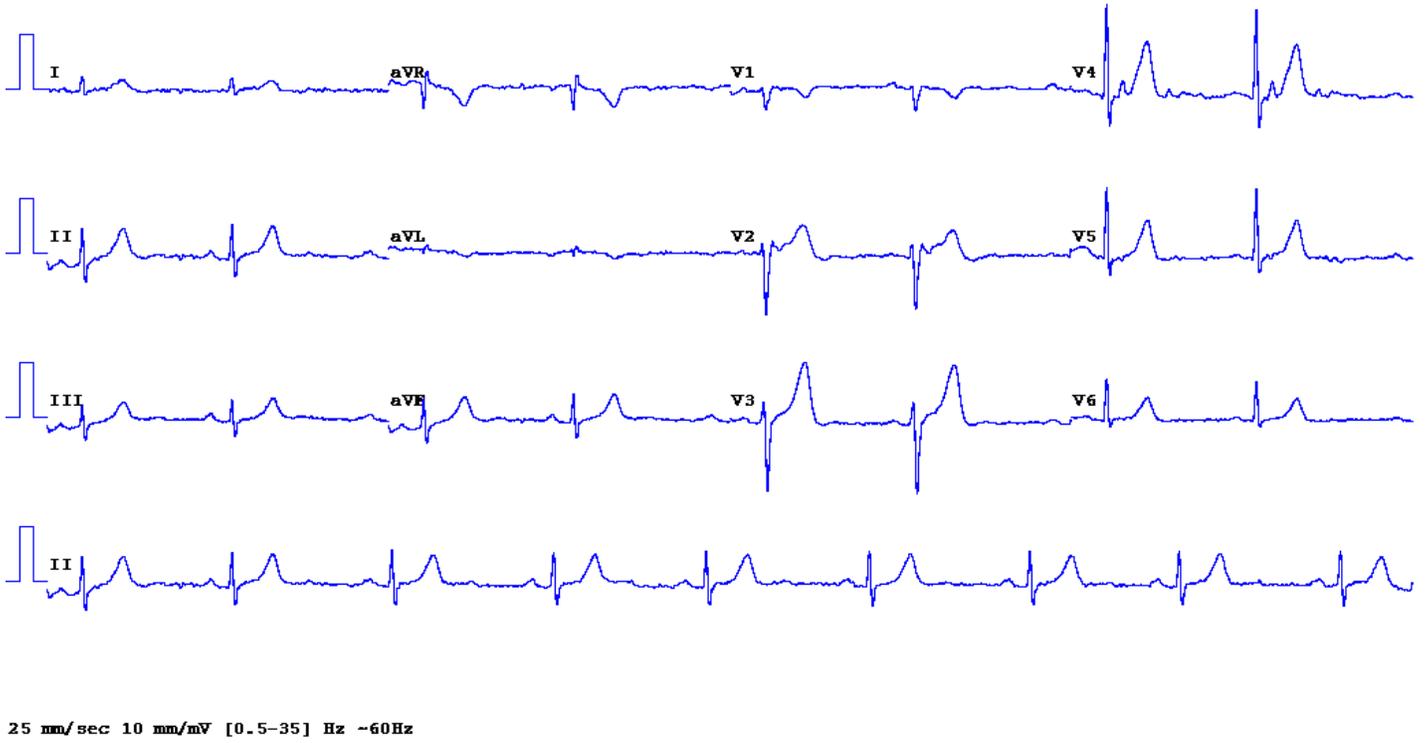

eFigure 17: Patient 09 LEAD 1 AVL STD ECG



eFigure 18: Patient 09 LEAD 1 AVL NEP ECG similar




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 19: Patient 10 LEAD 1 AVL STD ECG




25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 20: Patient 10 LEAD 1 AVL NEP ECG similar R waves

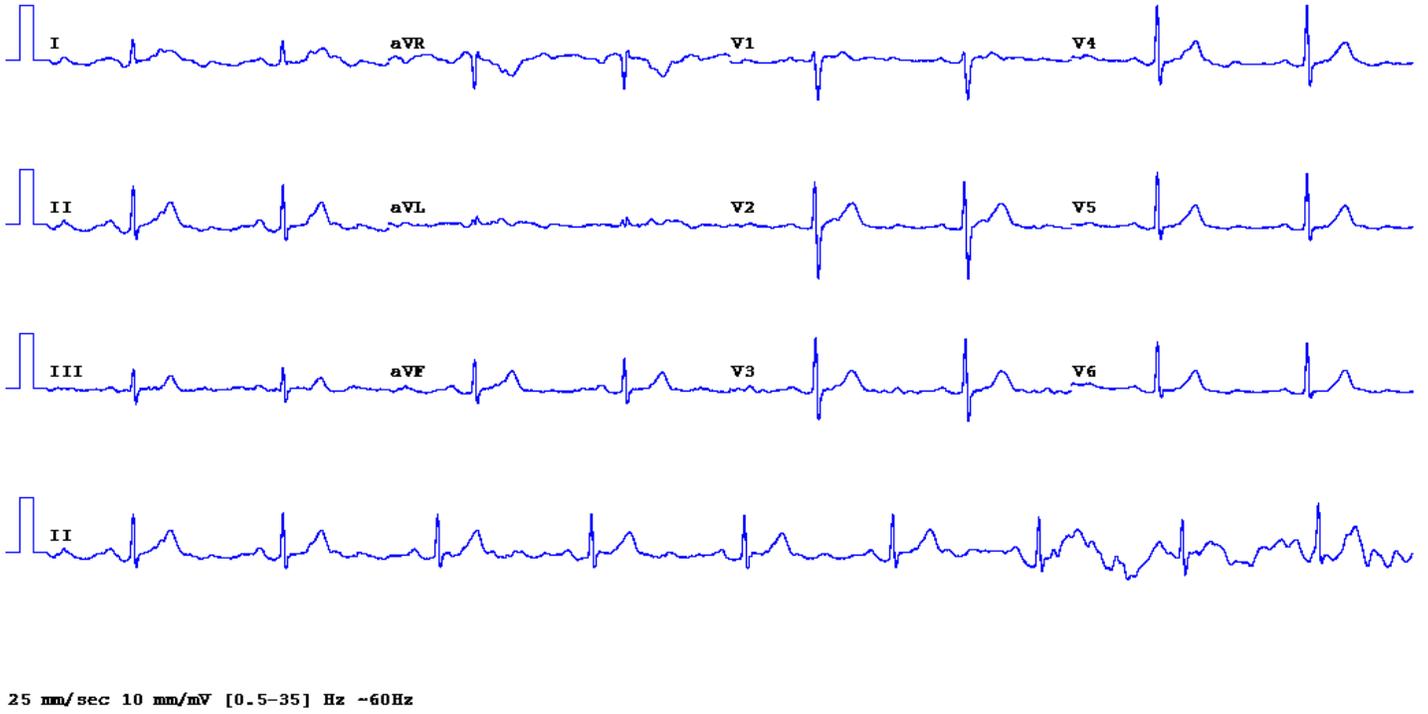

eFigure 21: Patient 11 LEAD 1 AVL STD ECG

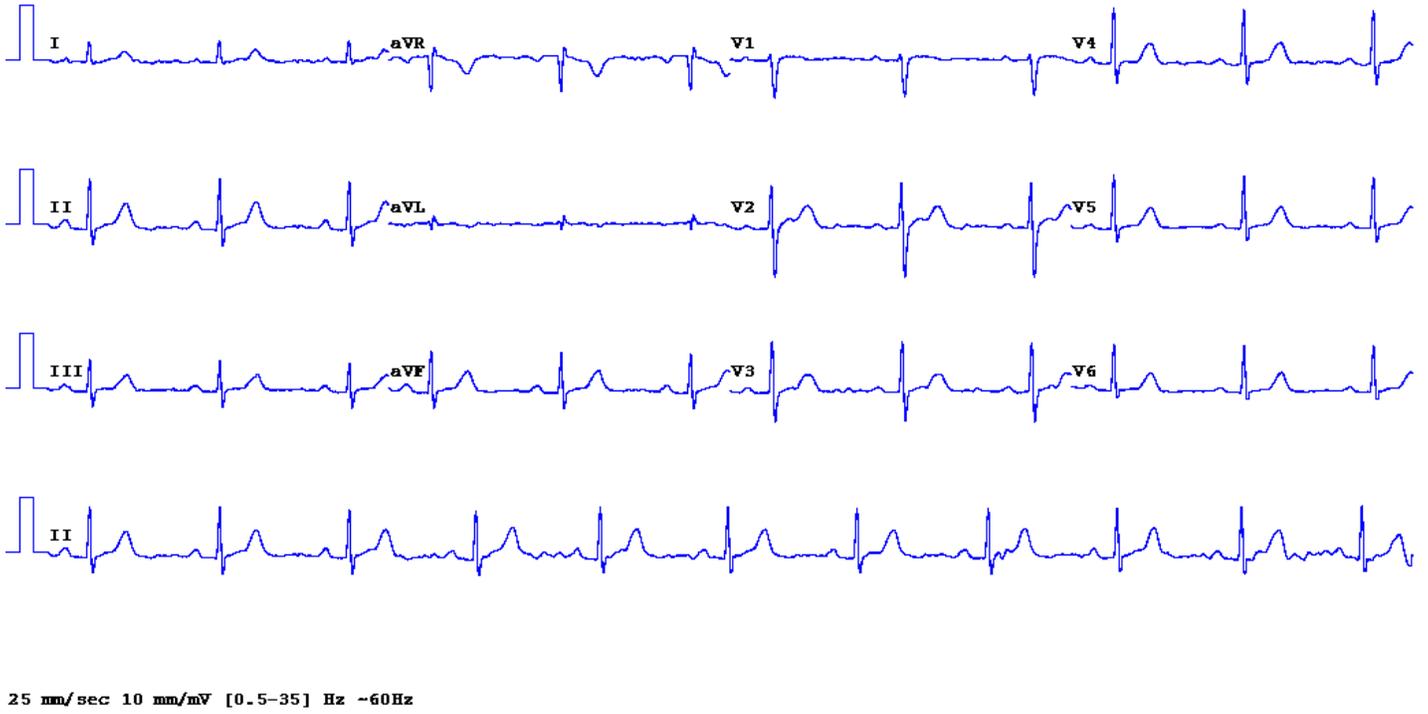

eFigure 22: Patient 11 LEAD 1 AVL NEP ECG all R waves amplitude similar to STD



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 23: Patient 12 LEAD 1 AVL STD ECG



25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 24: Patient 12 LEAD 1 AVL NEP ECG virtually identical.



eFigure 25: Patient 13 LEAD 1 AVL STD ECG




eFigure 26: Patient 13 LEAD 1 AVL NEP ECG R wave amplitude in inferior and other leads similar to STD



25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 27: Patient 14 LEAD 1 AVL STD ECG; very low R waves.

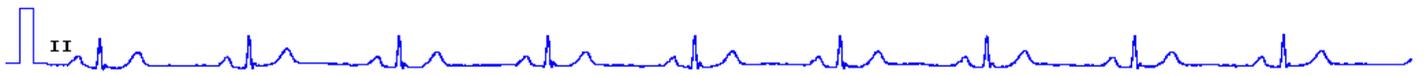

25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 28: Patient 14 LEAD 1 AVL NEP ECG similar; no QS emerged.




eFigure 29: Patient 15 LEAD 1 AVL STD ECG



eFigure 30: Patient 15 LEAD 1 AVL NEP ECG virtually identical R wave amplitude all leads.



eFigure 31: Patient 16 R lead 1 aVL STD ECG

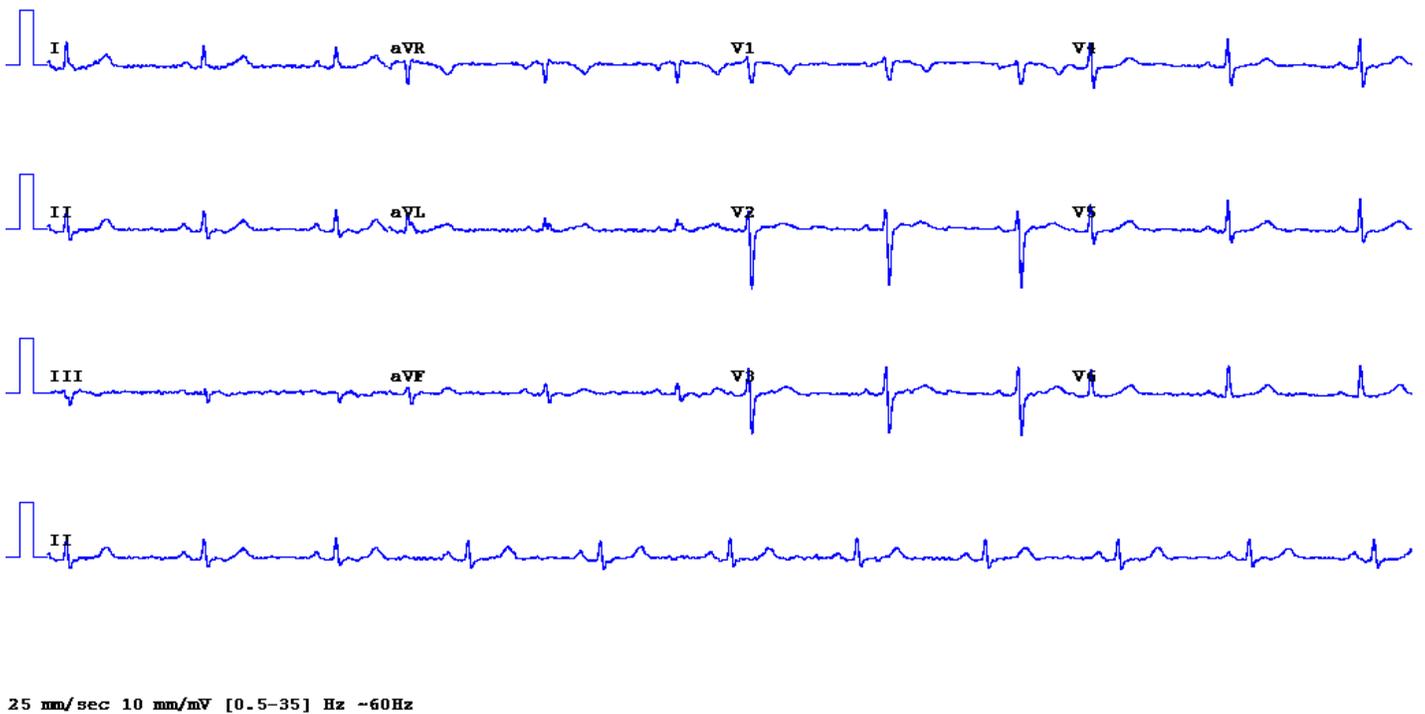

eFigure 32: Patient 16 R WAVE AMP NEP ECG NO INCREASE IN R WAVE AMPLITUDE IN 11 111 AVF AND NO DECREASE IN LEAD 1 AVL

## eAppendix 3




25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 1. 2014/05/30 Standard ECG, deep, wide inferior Q waves; old inferior infarct, QRS axis 19.




25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 2. NEP ECG shows virtually the same findings, deep, wide inferior Q waves: old inferior MI; QRS axis 25.

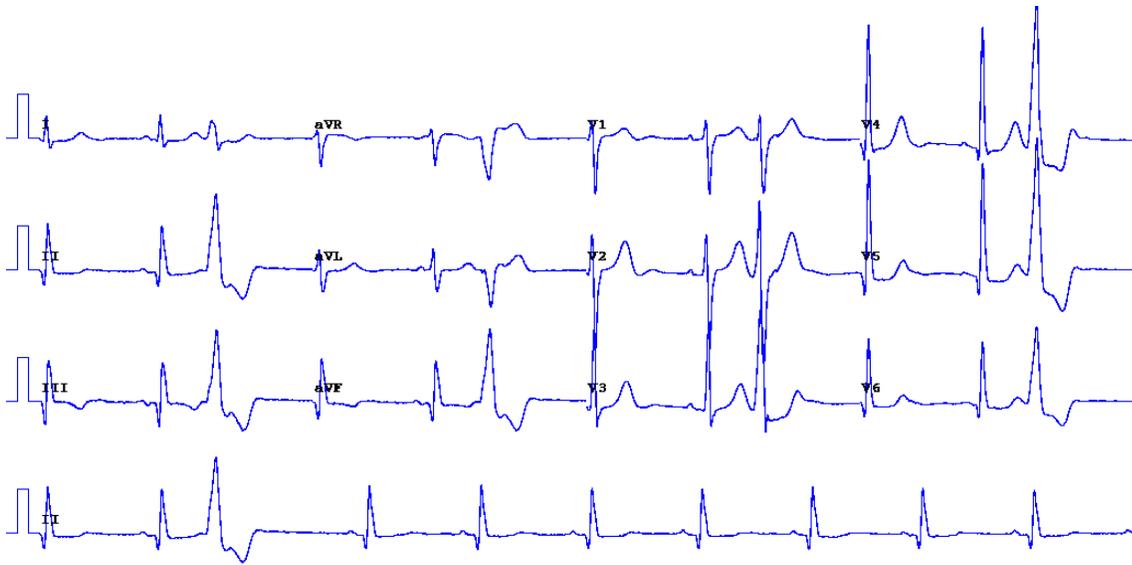

25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 3. 2014/ 05/09; standard ECG old inferior MI, QRS axis 74. VPBs.

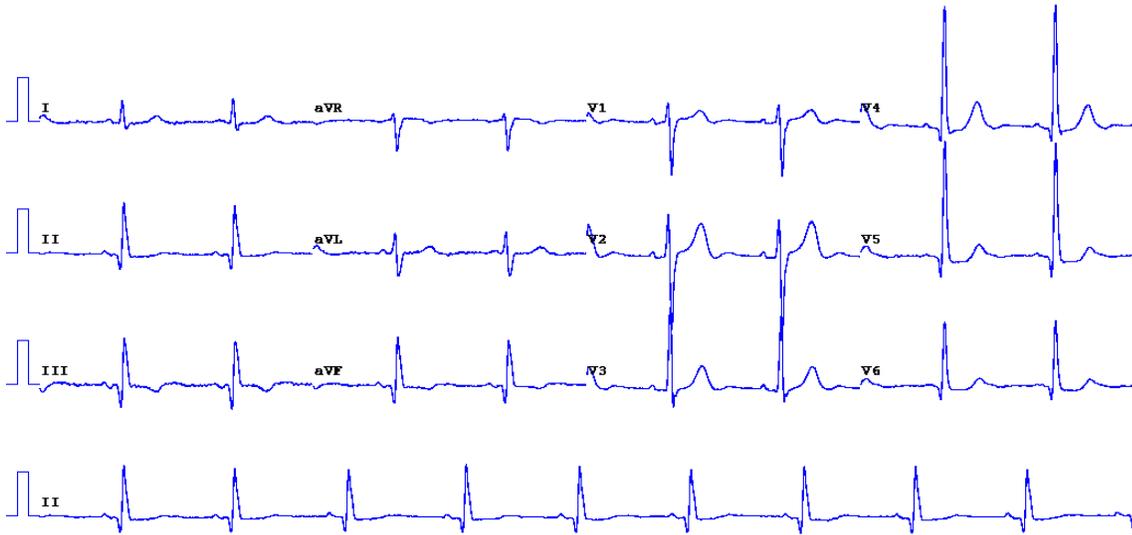

25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 4. new electrode placement ECG, similar Q waves in inferior leads, axis 69.



25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 5. standard ECG, 2014/05/23; Q waves 11,111, aVF: old inferior infarct, VPBs.



25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 6. NEP ECG reveals virtually the same, Q waves 11, 111, aVF, old inferior MI, VPBs.




25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 7. 2014 /05 /16 Standard 12 – lead ECG, Q waves 11, 111, aVF, old inferior MI



25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 8. New electrode placement ECG, similar wave form and inferior Q waves.

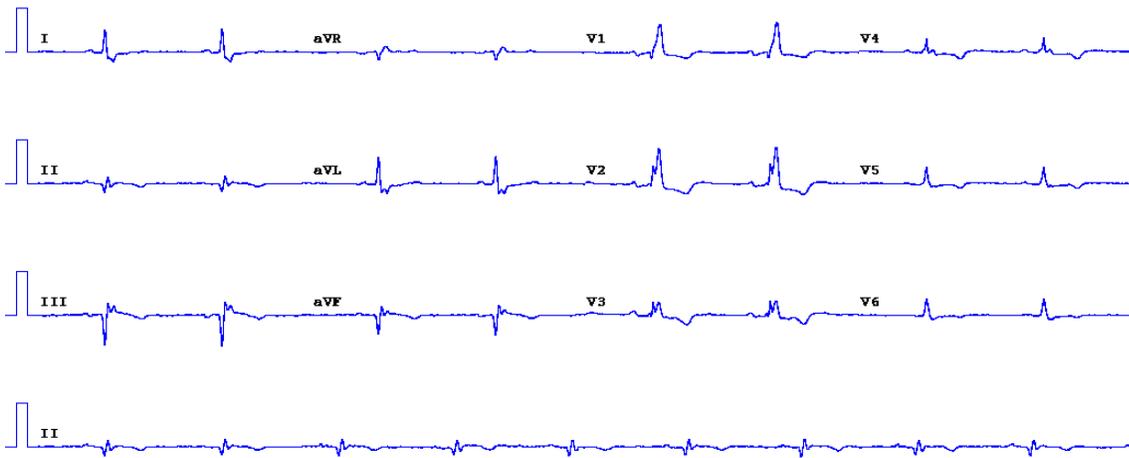


25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 9. 14/05/23 Standard ECG; old inferior MI, QRS axis  $-40$ . RBBB



25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 10. NEP ECG similar wave form, inferior Q waves, axis  $-35$ , RBBB. Virtually the same ST -T wave changes.

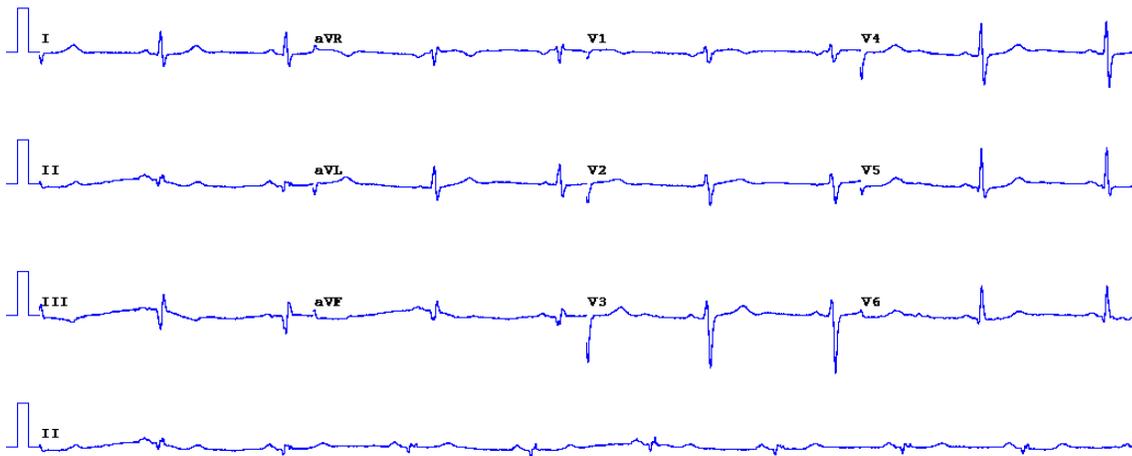

25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 11. 14/05/27; standard ECG; Q waves 11, 111 aVF; old inferior MI, axis -30



25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 12. NEP ECG similar inferior Q waves, axis -30; R wave amplitude in lead 1 aVL unchanged



eFigure 13. standard ECG, old inferior MI, axis -57 virtually the same wave form as NEP ECG shown below




eFigure 14. New electrode placement ECG; similar inferior Q waves; old inferior MI, axis -55



eFigure 15. standard ECG, old inferior MI; axis -47, note artifacts



eFigure 16. NEP ECG, old inferior MI similar to standard recording above, but artifacts are minimized and best appreciated if the tracing is enlarged; axis - 42

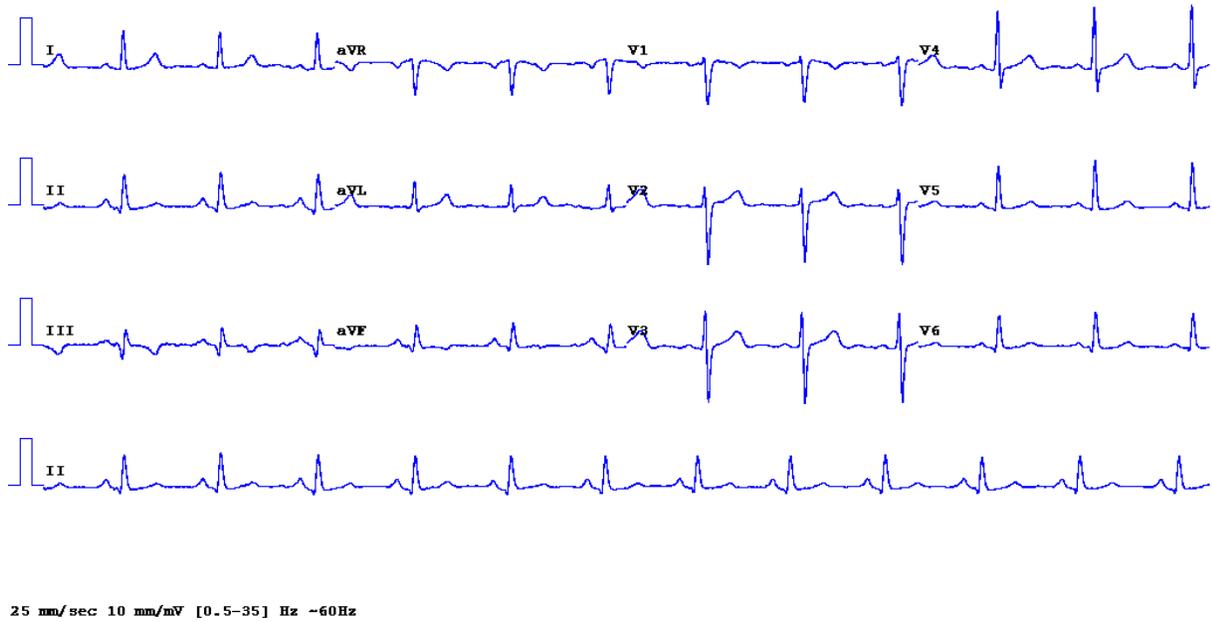

eFigure 17. Standard ECG, inferior Q waves: clinically proven old inferior MI



eFigure 18. NEP ECG, old inferior MI, virtually similar wave forms, inferior Q waves



25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 19. standard ECG; inferior Q waves; proven old inferior MI,




25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 20. NEP ECG, similar inferior Q waves; old inferior MI.




25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 21. Standard 12 – lead ECG; low amplitude R waves in leads 1 and aVL that can result in a QS pattern as observed with torso ( modified ) - leads

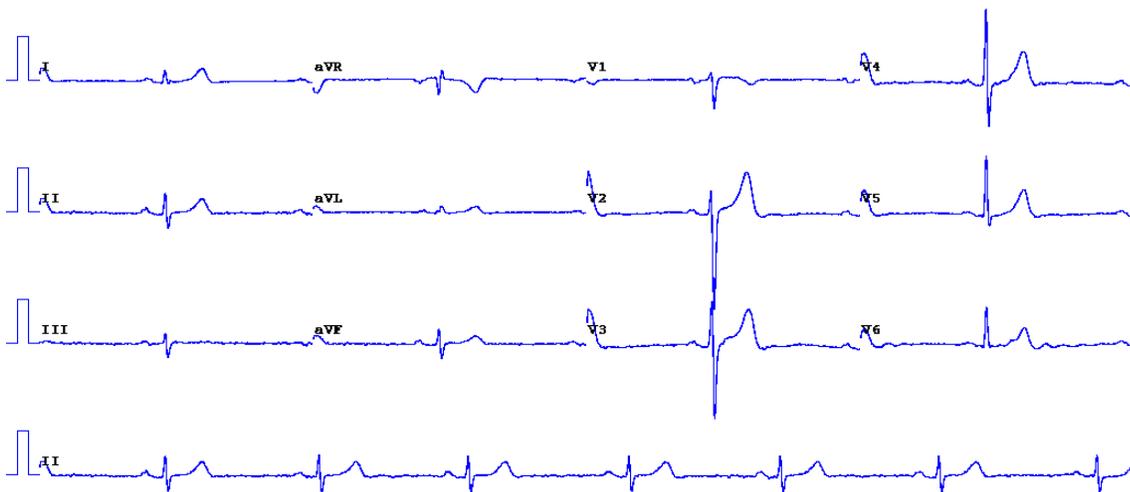


25 mm/sec 10 mm/mV [0.5-35] Hz -60Hz

eFigure 22. New electrode placement ECG showing no change in R wave amplitude in leads 1, aVL or inferior leads. Virtually the same wave forms.

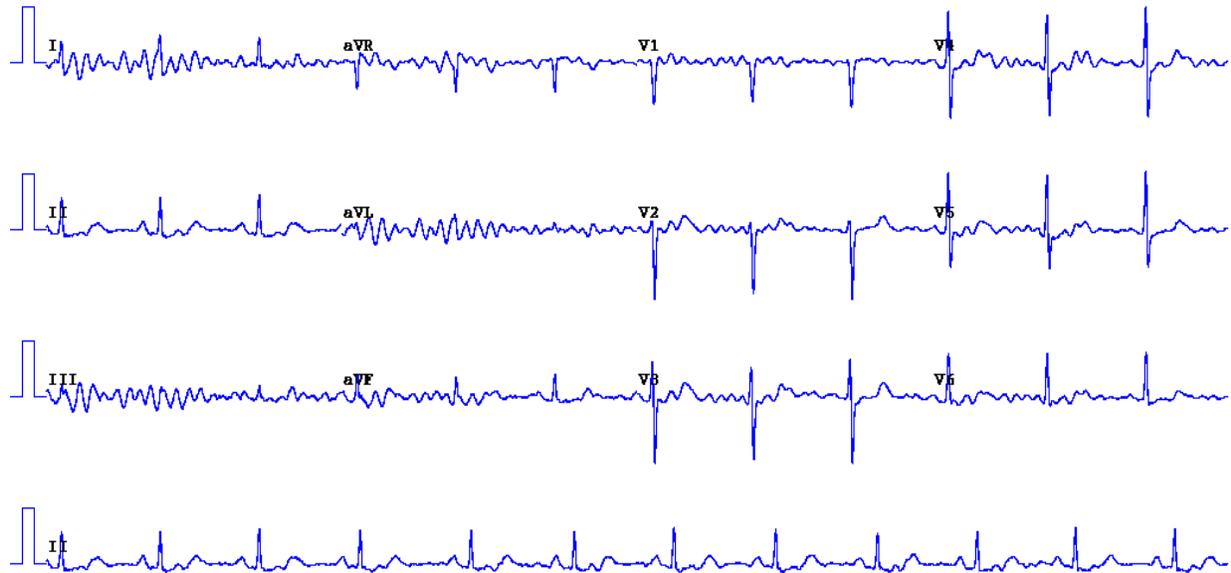

25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 23. July 2014. A standard method ECG shows artifacts; poor quality ECG. This problem we found to be more common in winter months; perhaps the wrist and ankle areas are colder than the torso and upper arm areas that allow better electrode sensing.

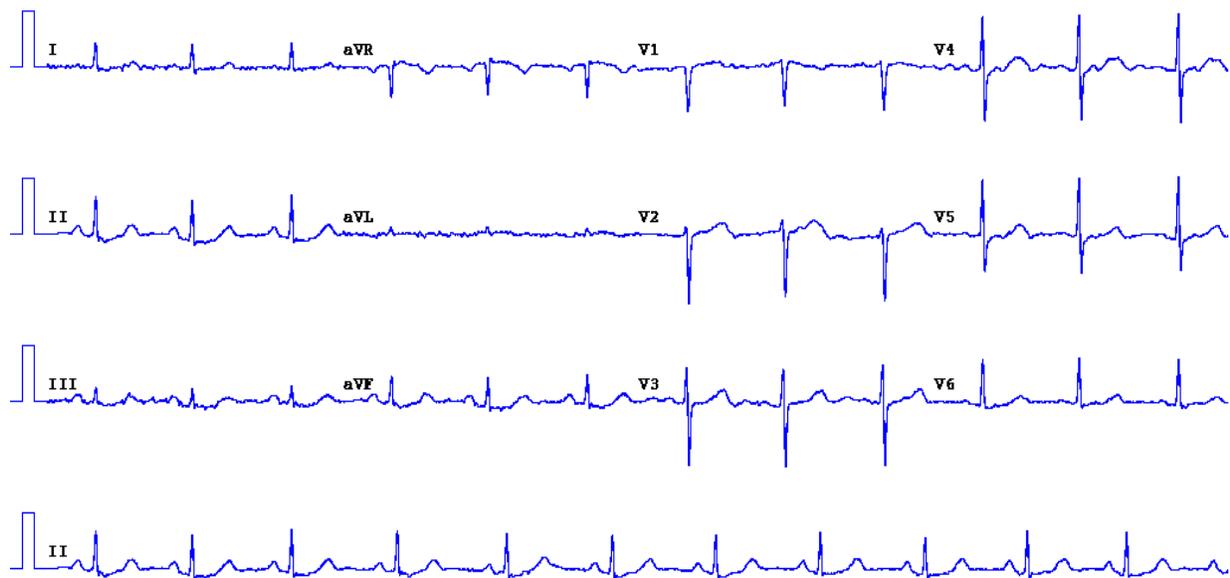

25 mm/sec 10 mm/mV [0.5-35] Hz ~60Hz

eFigure 24. NEP ECG done within 1 minute of STD shows clearing of artifacts